



Moving Optic Example:

- Velocity selector




Velocity Selectors


-

Select the neutron energy you want

VELOCITY SELECTORS


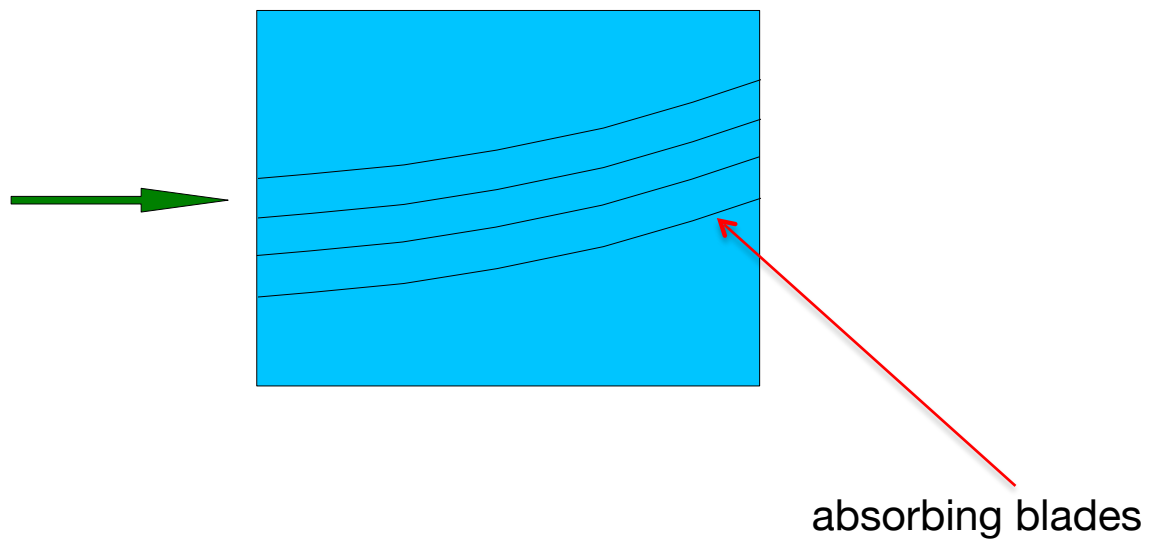


VELOCITY SELECTORS

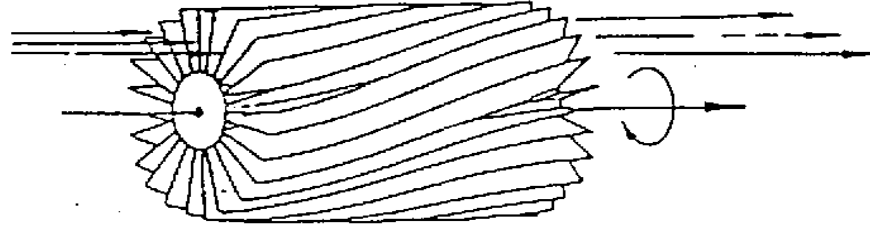


VELOCITY SELECTORS

‘broad’ monochromatization $\delta\lambda/\lambda \approx 10\%$

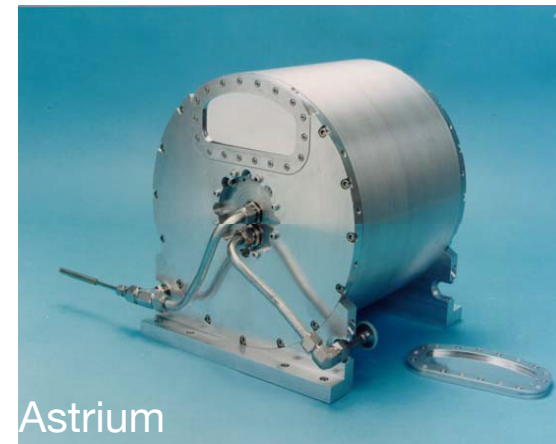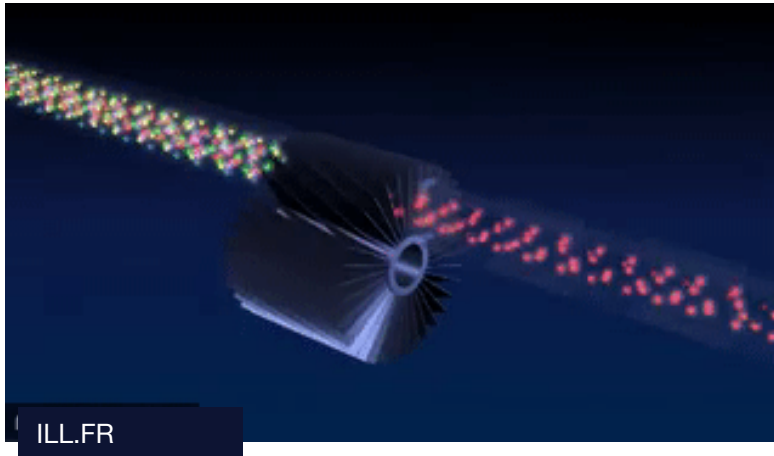

VELOCITY SELECTORS

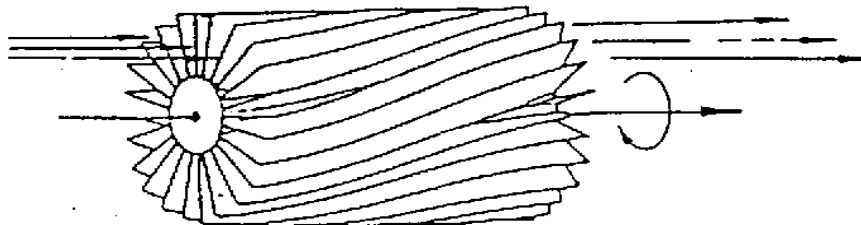


Example: `V_selector(xwidth=0.03, yheight=0.05, zdepth=0.30, radius=0.12, alpha=48.298, length=0.25, d=0.0004, nu=20000, nslit=72)`
 These are values for the D11@ILL Dornier 'Dolores' Velocity Selector (NVS 023)

Input parameters

Parameters in **boldface** are required; the others are optional.

Name	Unit	Description	Default
xwidth	m	Width of entry aperture	0.03
yheight	m	Height of entry aperture	0.05
zdepth	m	Distance between apertures, for housing containing the rotor	0.30
radius	m	Height from aperture centre to rotation axis	0.12
alpha	deg	Twist angle along the cylinder	48.298
length	m	Length of cylinder/rotor (less than zdepth)	0.25
d	m	Thickness of blades	0.0004
nu	Hz	Cylinder rotation speed, counter-clockwise, which is ideally $3956 \cdot \alpha \cdot \text{DEG2RAD} / 2 / \text{PI} / \lambda / \text{length}$	300
nslit	1	Number of Soller blades	72

VELOCITY SE

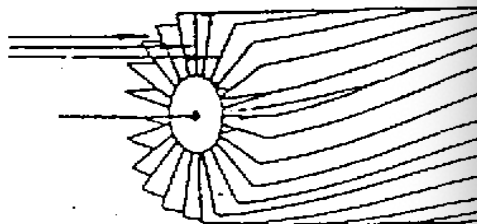

Example: `V_selector(xwidth=0.03, yheight=0.05, zdepth=0.30, radius=0.1)`
These are values for the D11@ILL Dornier 'Dolores' Velocity Selector

Input parameters

Parameters in **boldface** are required; the others are optional.

Name	Unit	Description
xwidth	m	Width of entry aperture
yheight	m	Height of entry aperture
zdepth	m	Distance between apertures, for housing containing the rot
radius	m	Height from aperture centre to rotation axis
alpha	deg	Twist angle along the cylinder
length	m	Length of cylinder/rotor (less than zdepth)
d	m	Thickness of blades
nu	Hz	Cylinder rotation speed, counter-clockwise, which is ideally $3956 \cdot \alpha \cdot \text{DEG2RAD} / 2 / \text{PI} / \lambda / \text{length}$
nslit	1	Number of Soller blades

