



Peter Willendrup

McStas samples for diffraction

Agenda

- General concepts (reminder)
- Focus on diffraction-oriented samples:
 - Incoherent scatterer
 - Powder
 - Single Crystal
 - File formats

Incoherent Sample

[Incoherent.comp](#)

Parameters in **boldface** are required; the others are optional.

An incoherent scatterer with various sample shape options

Name	Unit	Description	Default
geometry	str	Name of an Object File Format (OFF) or PLY file for complex geometry. The OFF/PLY file may be generated from XYZ coordinates using qhull/powercrust	0
radius	m	Outer radius of sample in (x,z) plane	0
xwidth	m	Horiz. dimension of sample (bounding box if off file), as a width	0
yheight	m	Vert. dimension of sample (bounding box if off file), as a height. A sphere shape is used when 0 and radius is set	0
zdepth	m	Depth of sample (bounding box if off file)	0
thickness	m	Thickness of hollow sample	0
target_x	-		0
target_y	m	position of target to focus at	0
target_z	-		0
focus_r	m	Radius of disk containing target. Use 0 for full space	0
focus_xw	m	horiz. dimension of a rectangular area	0
focus_yh	m	vert. dimension of a rectangular area	0
focus_aw	deg	horiz. angular dimension of a rectangular area	0
focus_ah	deg	vert. angular dimension of a rectangular area	0
target_index	1	Relative index of component to focus at, e.g. next is +1	0
pack	1	Packing factor	1
p_interact	1	MC Probability for scattering the ray; otherwise transmit	1
f_QE	1	Fraction of quasielastic scattering (rest is elastic)	0
gamma	1	Lorentzian width of quasielastic broadening (HWHM)	0
sigma_abs	barns	Absorption cross section pr. unit cell at 2200 m/s	5.08
sigma_inc	barns	Incoherent scattering cross section pr. unit cell	5.08
Vc	AA ³	Unit cell volume	13.827
concentric	1	Indicate that this component has a hollow geometry and may contain other components. It should then be duplicated after the inside part (only for box, cylinder, sphere)	0
order	-	Limit multiple scattering up to given order	0

Incoherent Sample

Incoherent.comp

Parameters in **boldface** are required; the others are optional.

Name	Unit	Description	Default
geometry	str	Name of an Object File Format (OFF) or PLY file for complex geometry. The OFF/PLY file may be generated from XYZ coordinates using qhull/powercrust	0
radius	m	Outer radius of sample in (x,z) plane	0
xwidth	m	Horiz. dimension of sample (bounding box if off file), as a width	0
yheight	m	Vert. dimension of sample (bounding box if off file), as a height. A sphere shape is used when 0 and radius is set	0
zdepth	m	Depth of sample (bounding box if off file)	0
thickness	m	Thickness of hollow sample	0
target_x	-		0
target_y	m	position of target to focus at	0
target_z	-		0
focus_r	m	Radius of disk containing target. Use 0 for full space	0
focus_xw	m	horiz. dimension of a rectangular area	0
focus_yh	m	vert. dimension of a rectangular area	
focus_aw	deg	horiz. angular dimension of a rectangular area	
focus_ah	deg	vert. angular dimension of a rectangular area	
target_index	1	Relative index of component in stack	
pack	1	Packing factor	
p_interact	1	MC Probability for scattering	
I_QE	1	Fraction of quasielastic intensity	
gamma	1	Lorentzian width of quasielastic peak	
sigma_abs	barns	Absorption cross section	
sigma_inc	barns	Incoherent scattering cross section	
Vc	AA ³	Unit cell volume	13.827
concentric	1	Indicate that this component has a hollow geometry and may contain other components. It should then be duplicated after the inside part (only for box, cylinder, sphere)	0
order	-	Limit multiple scattering up to given order	0

An incoherent scatterer with various sample shape options


p_interact = The probability that a given ray will interact (scatter) with a sample. Does not directly affect the intensity – only the quality of the resulting Monte Carlo estimate.

Crystalline Samples





From single crystal / crystallites to powder....


Single crystal

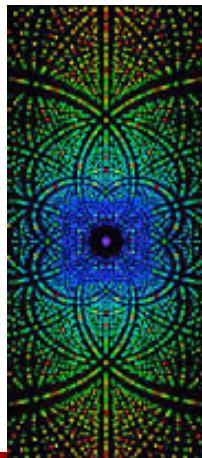


Polycrystal with a little disorder,
i.e. a *preferred orientation, texture*

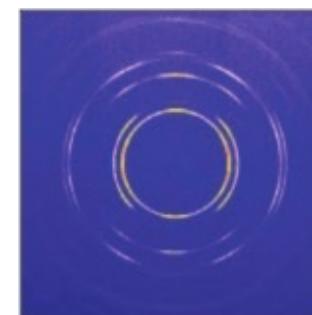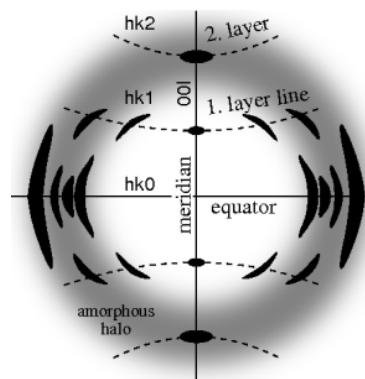

Powder with complete disorder

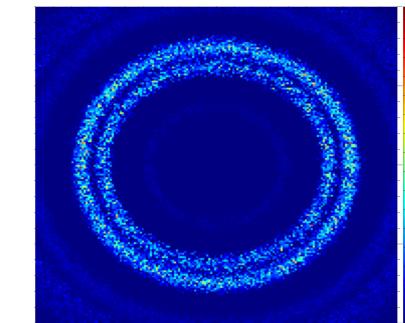

Single Crystal model

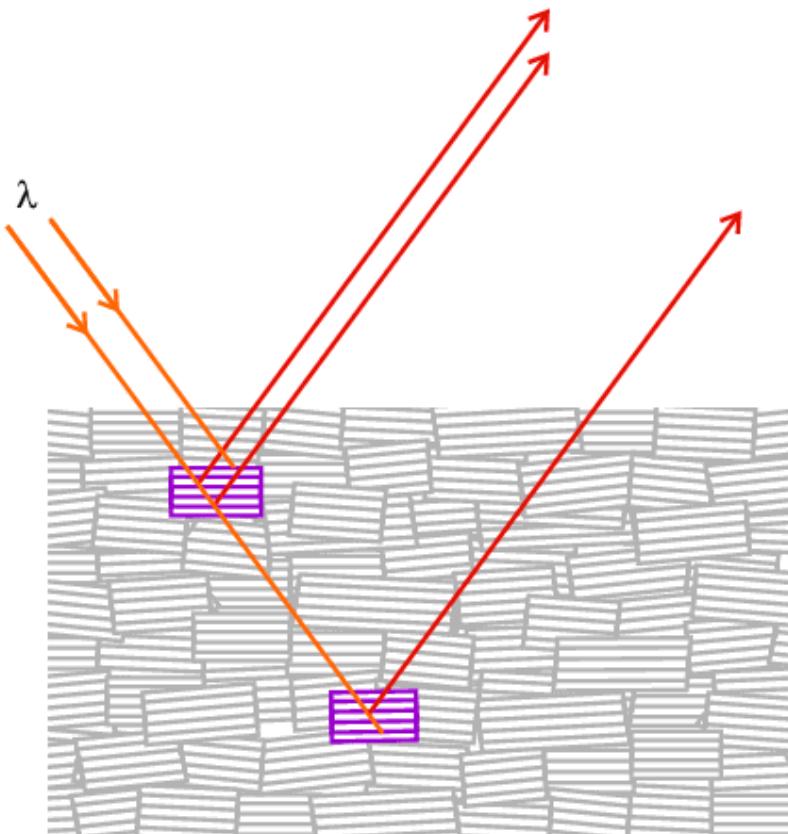

- *Models ideally imperfect crystal:*
- *Peak broadening is dominated by mosaic*
=> not appropriate for highly perfect crystals.
- “Space group symmetries” considered out of scope.
 - Optimized for small unit cells but...
 - No assumption about unit cell, orientation etc.
 - Gaussian model for mosaic (see later)

Courtesy University College London

Single Crystal

Input parameters

Parameters in **boldface** are required; the others are optional.

Name	Unit	Description	Default
mosaic_AB	arc_minutes, arc_minutes,1, 1, 1, 1, 1, 1	In Plane mosaic rotation and plane vectors (anisotropic), mosaic_A, mosaic_B, A_h,A_k,A_l, B_h,B_k,B_l. Puts the crystal in the in-plane mosaic state. Vectors A and B define plane in which the crystal roation is defined, and mosaic_A, mosaic_B, denotes the resp. mosaicities (gaussian RMS) with respect to the two reflections chosen by A and B (Miller indices).	Mosaic_AB_Undefined
reflections	string	File name containing structure factors of reflections. Use empty ("") or NULL for incoherent scattering only	0
geometry	str	Name of an Object File Format (OFF) or PLY file for complex geometry. The OFF/PLY file may be generated from XYZ coordinates using qhull/powercrust	0
xwidth	m	Width of crystal	0
yheight	m	Height of crystal	0
zdepth	no extinction simulated	[m] Depth of crystal	0
radius	m	Outer radius of sample in (x,z) plane	0
delta_d_d	1	Lattice spacing variance, gaussian RMS	1e-4
mosaic	arc minutes	Crystal mosaic (isotropic), gaussian RMS. Puts the crystal in the isotropic mosaic model state, thus disregarding other mosaicity parameters.	-1
mosaic_a	arc minutes	Horizontal (rotation around lattice vector a) mosaic (anisotropic), gaussian RMS. Put the crystal in the anisotropic crystal vector state. I.e. model mosaicity through rotation around the crystal lattice vectors. Has precedence over in-plane mosaic model.	-1
mosaic_b	arc minutes	Vertical (rotation around lattice vector b) mosaic (anisotropic), gaussian RMS.	-1
mosaic_c	arc minutes	Out-of-plane (Rotation around lattice vector c) mosaic (anisotropic), gaussian RMS	-1
recip_cell	1	Choice of direct/reciprocal (0/1) unit cell definition	0
barns	1	Flag to indicate if $ F ^2$ from 'reflections' is in barns or fm ² . barns=1 for laz and isotropic constant elastic scattering (reflections=NULL), barns=0 for lau type files	0
ax	AA or AA ⁻¹	Coordinates of first (direct/recip) unit cell vector	0
ay			0
az			0
bx	AA or AA ⁻¹	Coordinates of second (direct/recip) unit cell vector	0
by			0
bz			0
cx	AA or AA ⁻¹	Coordinates of third (direct/recip) unit cell vector	0
cy			0
cz			0
p_transmit	1	Monte Carlo probability for neutrons to be transmitted without any scattering. Used to improve statistics from weak reflections	-1
sigma_abs	barns	Absorption cross-section per unit cell at 2200 m/s	0
sigma_inc	barns	Incoherent scattering cross-section per unit cell Use -1 to unactivate	0
aa	deg	Unit cell angles alpha, beta and gamma. Then uses norms of vectors a,b and c as lattice parameters	0
bb	deg	Beta angle	0
cc	deg	Gamma angle	0
order	0: all, 1: first, 2: second, ...	[1] Limit multiple scattering up to given order	0
RX	m	Radius of horizontal along X lattice curvature. flat for 0	0
RY	m	Radius of vertical lattice curvature. flat for 0	0
RZ	m	Radius of horizontal along Z lattice curvature. flat for 0	0
powder	1	Flag to indicate powder mode, for simulation of Debye-Scherrer cones via random crystallite orientation. A powder texture can be approximated with 0	0
PG	1	Flag to indicate "Pyrolytic Graphite" mode, only meaningful with choice of Graphite.lau, models PG crystal. A powder texture can be approximated with 0	0

Single Crystal

Input parameters

Parameters in **boldface** are required; the others are optional.

Name	Unit	Description	Default
mosaic_AB	arc_minutes, arc_minutes, 1, 1, 1, 1, 1	In Plane mosaic rotation and plane vectors (anisotropic), mosaic_A, mosaic_B, A_h,A_k,A_l, B_h,B_k,B_l. Puts the crystal in the in-plane mosaic state. Vectors A and B define plane in which the crystal roation is defined, and mosaic_A, mosaic_B, denotes the resp. mosaicities (gaussian RMS) with respect to the two reflections chosen by A and B (Miller indices).	Mosaic_AB_Undefined
reflections	string	File name containing structure factors of reflections. Use empty ("") or NULL for incoherent scattering only	0
geometry	str	Name of an Object File Format (OFF) or PLY file for complex geometry. The OFF/PLY file may be generated from XYZ coordinates using qhull/powercrust	0
xwidth	m	Width of crystal	0
yheight	m	Height of crystal	0
zdepth	no extinction simulated	[m] Depth of crystal	0
radius	m	Outer radius of sample in (x,z) plane	0
delta_d_d	1	Lattice spacing variance, gaussian RMS	1e-4
mosaic	arc minutes	Crystal mosaic (isotropic), gaussian RMS. Puts the crystal in the isotropic mosaic model state, thus disregarding other mosaicity parameters.	-1
mosaic_a	arc minutes	Horizontal (rotation around lattice vector a) mosaic (anisotropic), gaussian RMS. Put the crystal in the anisotropic crystal vector state. I.e. model mosaicity through rotation around the crystal lattice vectors. Has precedence over in-plane mosaic model.	-1
mosaic_b	arc minutes	Vertical (rotation around lattice vector b) mosaic (anisotropic), gaussian RMS.	-1
mosaic_c	arc minutes	Out-of-plane (Rotation around lattice vector c) mosaic (anisotropic), gaussian RMS	-1
recip_cell	1	Choice of direct/reciprocal (0/1) unit cell definition	0
barns	1	Flag to indicate if $ F ^2$ from 'reflections' is in barns or fm^2. barns=1 for laz and isotropic constant elastic scattering (reflections=NULL), barns=0 for lau type files	0
ax	AA or AA^-1	Coordinates of first (direct/recip) unit cell vector	0
ay			0
az			0
bx	AA or AA^-1	Coordinates of second (direct/recip) unit cell vector	0
by			0
bz			0
cx	AA or AA^-1		0
cy			0
cz			0
p_transm	1		-1
sigma_pos	barns		0
sigm_a_inc	barns	Incoherent scattering cross-section per unit cell Use -1 to deactivate	0
aa	deg	Unit cell angles alpha, beta and gamma. Then uses norms of vectors a,b and c as lattice parameters	0
bb	deg	Beta angle	0
cc	deg	Gamma angle	0
order	0: all, 1: first, 2: second, ...	[1] Limit multiple scattering up to given order	0
RX	m	Radius of horizontal along X lattice curvature. flat for 0	0
RY	m	Radius of vertical lattice curvature. flat for 0	0
RZ	m	Radius of horizontal along Z lattice curvature. flat for 0	0
powder	1	Flag to indicate powder mode, for simulation of Debye-Scherrer cones via random crystallite orientation. A powder texture can be approximated with 0	0
PG	1	Flag to indicate "Pyrolytic Graphite" mode, only meaningful with choice of Graphite.lau, models PG crystal. A powder texture can be approximated with 0	0

Unit cell vectors (ax,ay,az, bx, by,bz, cx,cy, cz). Encodes the crystal orientation relative to the Component's orientation

Single Crystal


Input parameters

Parameters in **boldface** are required; the others are optional.


Name	Unit	Description	Default
mosaic_AB	arc_minutes, arc_minutes,1, 1, 1	In Plane mosaic rotation and plane vectors (anisotropic), mosaic_A, mosaic_B, A_h.A_k,A_l, B_h.B_k,B_l. Puts the crystal in the in-plane mosaic state. Vectors A and B define plane in which the crystal roation is defined, and mosaic_A, mosaic_B, denotes the resp. mosaicities (gaussian RMS) with respect to the two reflections chosen by A and B (Miller indices).	Mosaic_AB_Undefined
reflections	string	File name containing structure factors of reflections. Use empty ("") or NULL for incoherent scattering only	0
g_mosaic	str	Name of an Object File Format (OFF) or PLY file for complex geometry. The OFF/PLY file may be generated from XYZ coordinates using qhull/powercrust	0
xwidth	m	Width of crystal	0
yheight	m	Height of crystal	0
zdepth	no extinction simulated	[m] Depth of crystal	0
radius	m	Outer radius of sample in (x,z) plane	0
lens_d	1		1e-4
mosaic	arc minutes		-1
mosaic_a	arc minutes		-1
mosaic_b	arc minutes		-1
mosaic_c	arc minutes		-1
recip_cell	1		0
barns	1	Flag to indicate if $ F ^2$ from 'reflections' is in barns or fm 2 . barns=1 for laz and isotropic constant elastic scattering (reflections=NULL), barns=0 for lau type files	0
cx	AA or AA $^{-1}$	Coordinates of first (direct/recip) unit cell vector	0
ay			0
az			0
bx	AA or AA $^{-1}$	Coordinates of second (direct/recip) unit cell vector	0
by			0
bz			0
cx	AA or AA $^{-1}$	Coordinates of third (direct/recip) unit cell vector	0
cy			0
cz			0
p_transmit	1	Monte Carlo probability for neutrons to be transmitted without any scattering. Used to improve statistics from weak reflections	-1
sigma_abs	barns	Absorption cross-section per unit cell at 2200 m/s	0
sigma_inc	barns	Incoherent scattering cross-section per unit cell Use -1 to unactivate	0
aa	deg	Unit cell angles alpha, beta and gamma. Then uses norms of vectors a,b and c as lattice parameters	0
bb	deg	Beta angle	0
cc	deg	Gamma angle	0
order	0, all, 1: first, 2: second, ...	[1] Limit multiple scattering up to given order	0
RX	m	Radius of horizontal along X lattice curvature. flat for 0	0
RY	m	Radius of vertical lattice curvature. flat for 0	0
RZ	m	Radius of horizontal along Z lattice curvature. flat for 0	0
powder	1	Flag to indicate powder mode, for simulation of Debye-Scherrer cones via random crystallite orientation. A powder texture can be approximated with 0	0
PG	1	Flag to indicate "Pyrolytic Graphite" mode, only meaningful with choice of Graphite.lau, models PG crystal. A powder texture can be approximated with 0	0

Mosaicity parameters: mosaic, mosaic_a, mosaic_b, mosaic_c, mosaic_AB

Single Crystal mosaicity

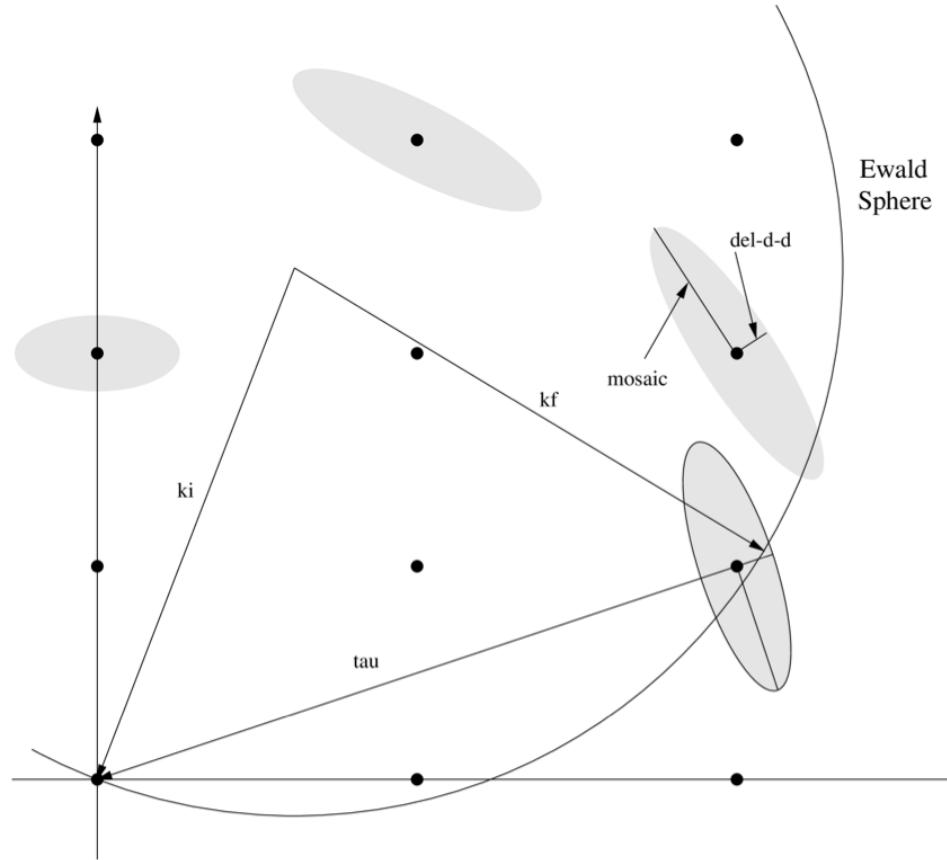


Single Crystal mosaicity

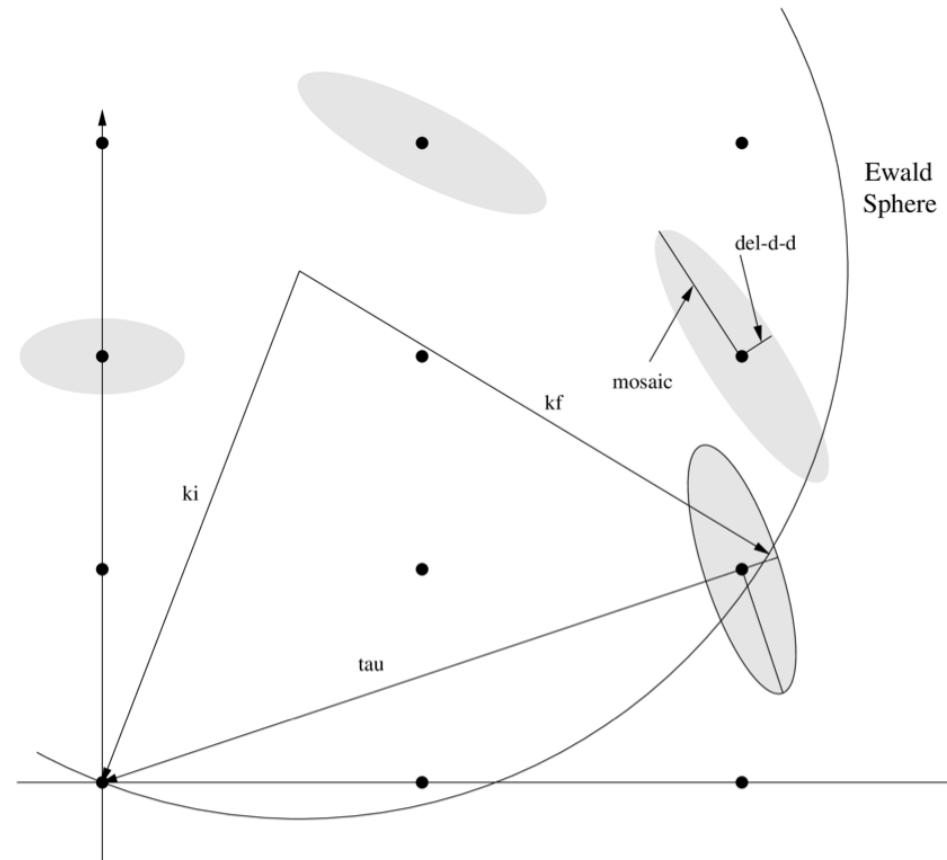

For an arbitrary, incoming particle:

- * Lay out the Ewald construction
(assume elastic scattering)
- * Find candidate reflections
(widened from Dirac-points in rcp. space by mosaicity and lattice space variation $\delta\mathbf{d} / \mathbf{d}$)
- * Perform Monte Carlo choices and calculate reflection probability, taking all coherent, incoherent, absorption cross-sections into account
- * Define an outgoing particle

Single Crystal mosaicity

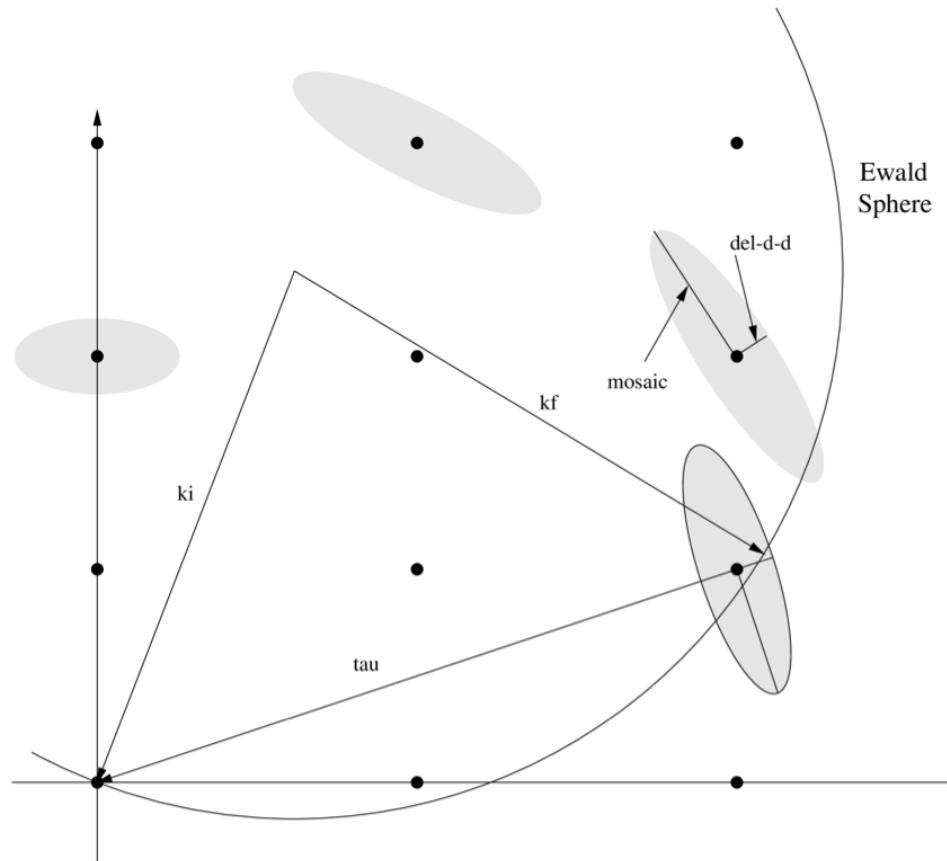

For an arbitrary, incoming particle:


- * Lay out the Ewald construction
(assume elastic scattering)
- * Find candidate reflections
(widened from Dirac-points in rcp. space by mosaicity and lattice space variation $\delta d / d$)
- * Perform Monte Carlo choices and calculate reflection probability, taking all coherent, incoherent, absorption cross-sections into account
- * Define an outgoing particle

Repeat until exiting crystal (*)

Single Crystal mosaicity

Three models:



- Isotropic Mosaicity = Gaussian distribution of rotation angle, with given standard deviation, around a random vector.
- Anistropic mosaicity = Gaussian distribution of angles around each of the unit cell axis.

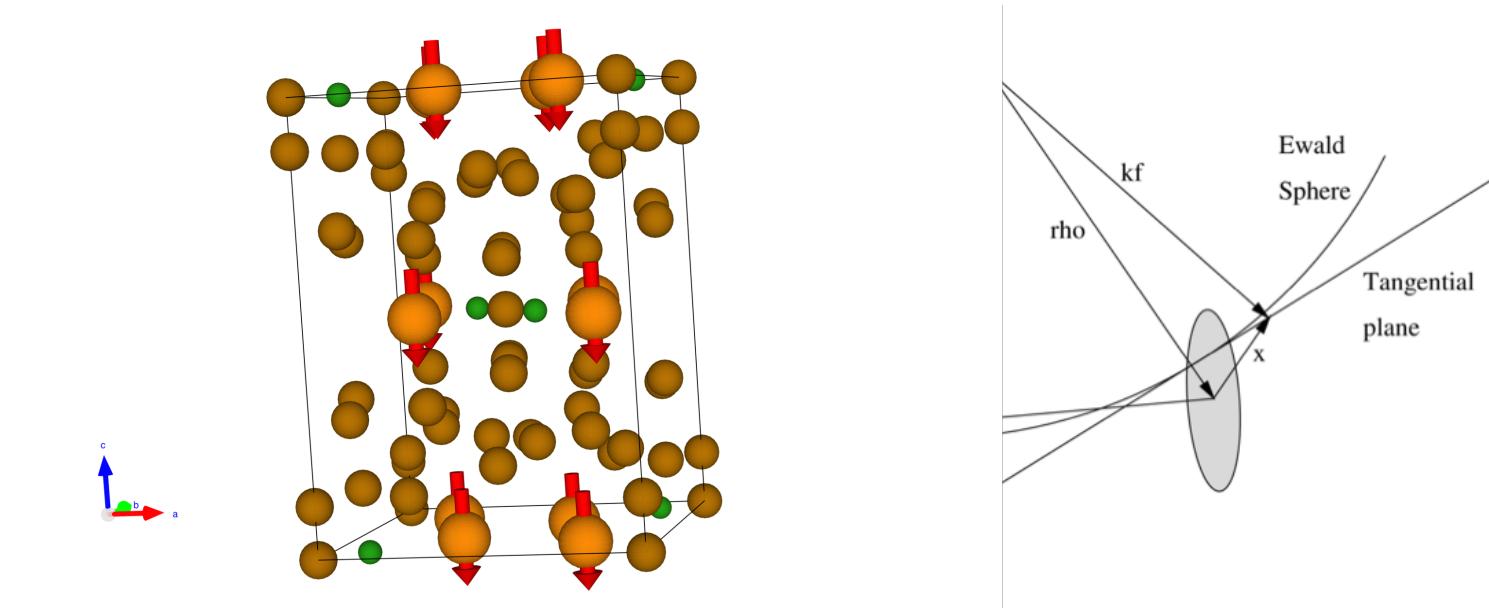

Single crystal mosaicity

Three models:

- Isotropic Mosaicity = Gaussian distribution of rotation angle, with given standard deviation, around a random vector.
- Anistropic mosaicity = Gaussian distribution of angles around each of the unit cell axis.
- “Rocking curve mosaic” = Gaussian distr. around axis determined by width of two given Bragg peaks. Specify σ , and hkl for two peaks.

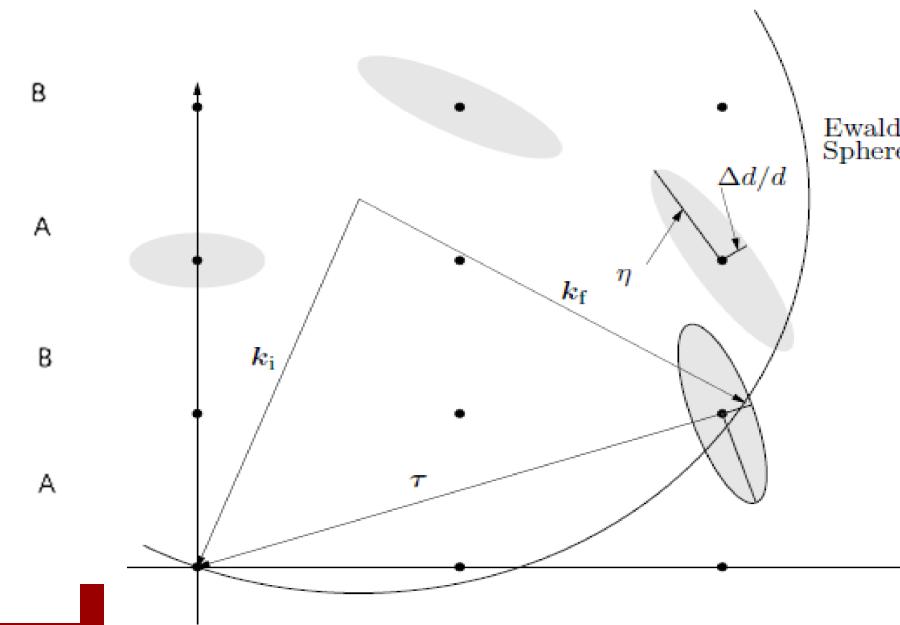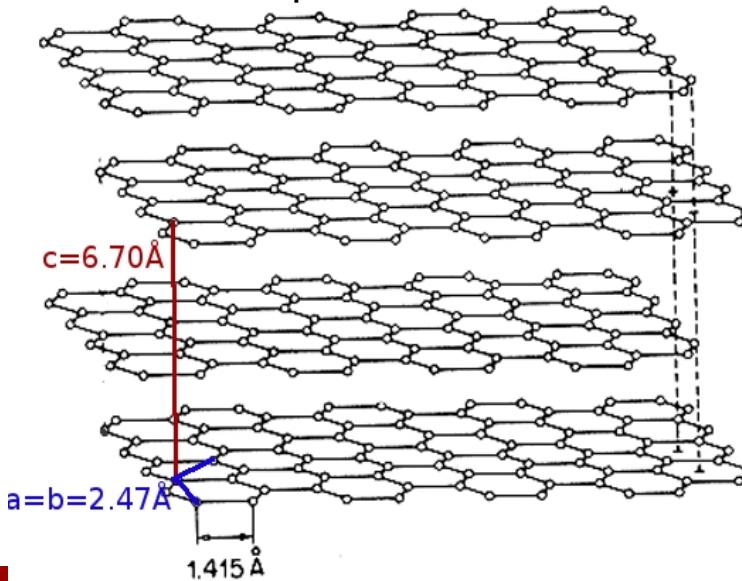

9.2 PowderN

Among the most common materials measured with neutron scattering techniques are powders, which are composed of a large number of tiny single crystals. The scattering intensity, by means of constructive interference and averaging over all crystallites, appears as well defined rings around the out-going direct beam direction. These rings are also present in liquids, but are much smoother, as we shall see in next section.

Let us consider a material of density ρ made of N scattering units each with a unit cell volume V_0 arranged with lattice spacings $d_Q=2\pi/Q$ associated to structure factors $F(Q)$ with equivalent reflection multiplicities j_Q . These structure factors characterize the efficiency of the reflection with momentum exchange Q . Following Squires [14] the scattering probability for an incoming neutron with wavelength λ penetrating along a distance x into the material is about $1-\exp(-\rho\sigma_{cone}x)$ where σ_{cone} is the so-called coherent elastic cross-section of the ring,

$$\sigma_{cone} = \frac{N \pi \lambda^2 j_Q |F(Q)|^2}{V_0 Q} .$$

This relation is only valid under certain conditions, among which $d_Q > \lambda/2$, from the Bragg law. As the possible lattice spacings d in the material can not exceed a maximum value (for instance the inter-atomic distance), it appears that as the neutron wavelength increases, the number of visible rings in the diffractogram will decrease, until no more scattering is possible above the so-called *Bragg edge*, where materials become transparent to neutrons (except for absorption and incoherent scattering). This is why most of the diffractometers use thermal and hot neutrons. Cold neutrons can only scatter on large distance arrangements in materials, *e.g.* in larger molecules and proteins.

PowderN inputs

Input parameters

Parameters in **boldface** are required; the others are optional.

Name	Unit	Description	Default
format	no quotes	Name of the format, or list of column indexes (see Description).	Undefined
reflections			"NULL"
geometry	str	Name of an Object File Format (OFF) or PLY file for complex geometry. The OFF/PLY file may be generated from XYZ coordinates using qhull/powercrust	"NULL"
radius	m	Outer radius of sample in (x,z) plane	0
yheight	m	Height of sample y direction	0
xwidth	m	Horiz. dimension of sample, as a width	0
zdepth	m	Depth of box sample	0
thickness			0
pack	1	Packing factor	1
Vc	AA^3	Volume of unit cell=nb atoms per cell/density of atoms	0
sigma_abs	barns	Absorption cross section per unit cell at 2200 m/s. Use a negative value to unactivate it	0
sigma_inc	barns	Incoherent cross section per unit cell. Use a negative value to unactivate it	0
delta_d_d	0/1	Global relative delta_d_d/d broadening when the 'w' column is not available. Use 0 if ideal.	0
p_inc	1	Fraction of incoherently scattered neutron rays	0.1
p_transmit	1	Fraction of transmitted (only attenuated) neutron rays	0.1
DW	1	Global Debye-Waller factor when the 'DW' column is not available. Use 1 if included in F2	0
nb_atoms	1	Number of sub-unit per unit cell, that is ratio of sigma for chemical formula to sigma per unit cell	1
d_phi	deg	Angle corresponding to the vertical angular range to focus to, e.g. detector height. 0 for no focusing	0
p_interact	1	Fraction of events interacting with sample, e.g. 1-p_transmit-p_inc	0
concentric	only for box, cylinder, sphere	[1] Indicate that this component has a hollow geometry and may contain other components. It should then be duplicated after the inside part	0
density	g/cm^3	Density of material. rho=density/weight/1e24*N_A.	0
weight	g/mol	Atomic/molecular weight of material	0
barns	1	Flag to indicate if F ^2 from 'reflections' is in barns or fm^2 (barns=1 for laz, barns=0 for lau type files).	1
Strain	ppm	Global relative delta_d_d/d shift when the 'Strain' column is not available. Use 0 if ideal.	0
focus_flip	1	Controls the sense of d_phi. If 0 d_phi is measured against the xz-plane. If !=0 d_phi is measured against zy-plane.	0

PowderN inputs

Input parameters

Parameters in **boldface** are required; the others are optional.

Name	Unit	Description	Default
format	no quotes	Name of the format, or list of column indexes (see Description).	Undefined
reflections			"NULL"
geometry	str	Name of an Object File Format (OFF) or PLY file for complex geometry. The OFF/PLY file may be generated from XYZ coordinates using qhull/powercrust	"NULL"
radius	m	Outer radius of sample in (x,z) plane	0
yheight	m	Height of sample y direction	0
xwidth	m	Horiz. dimension of sample, as a width	0
zdepth	m	Depth of box sample	0
thickness			0
pack	1	Packing factor	1
Vc	AA^3	Volume of unit cell=nb atoms per cell/density of atoms	0
sigma_abs	barns	Absorption cross section per unit cell at 2200 m/s. Use a negative value to unactivate it	0
sigma_inc	barns	Incoherent cross section per unit cell. Use a negative value to unactivate it	0
delta_d_d	0/1		0
p_inc	1		0.1
p_transmit	1	Fraction of transmitted (only attenuated) neutron rays	0.1
DW	1	Global Debye-Waller factor when the 'DW' column is not available. Use 1 if included in F2	0
nb_atoms	1	Number of sub-unit per unit cell, that is ratio of sigma for chemical formula to sigma per unit cell	1
d_phi	deg	Angle corresponding to the vertical angular range to focus to, e.g. detector height. 0 for no focusing	0
p_interact	1	Fraction of events interacting with sample, e.g. 1-p_transmit-p_inc	0
concentric	only for box, cylinder, sphere	[1] Indicate that this component has a hollow geometry and may contain other components. It should then be duplicated after the inside part	0
density	g/cm^3	Density of material. rho=density/weight/1e24*N_A.	0
weight	g/mol	Atomic/molecular weight of material	0
barns	1	Flag to indicate if F ^2 from 'reflections' is in barns or fm^2 (barns=1 for laz, barns=0 for lau type files).	1
Strain	ppm	Global relative delta_d_d shift when the 'Strain' column is not available. Use 0 if ideal.	0
focus_flip	1	Controls the sense of d_phi. If 0 d_phi is measured against the xz-plane. If !=0 d_phi is measured against zy-plane.	0

Sampling parameters

Concentric = Hollow
Powder [cylinder]


PowderN inputs

Input parameters

Parameters in **boldface** are required; the others are optional.

Name	Unit	Description
format	no quotes	Name of the format, or list of column indexes (see Description).
reflections		
geometry	str	Name of an Object File Format (OFF) or PLY file for complex geometry. The OFF/PLY file may be generated from XYZ coordinates us
radius	m	Outer radius of sample in (x,z) plane
yheight	m	Height of sample y direction
xwidth	m	Horiz. dimension of sample, as a width
zdepth	m	Depth of box sample
thickness		
pack	1	Packing factor
Vc	AA^3	Volume of unit cell=nb atoms per cell/density of atoms
sigma_abs	barns	Absorption cross section per unit cell at 2200 m/s. Use a negative value to unactivate it
sigma_inc	barns	Incoherent cross section per unit cell. Use a negative value to unactivate it
delta_d_d	0/1	In the 'w' column is not available. Use 0 if ideal.
p_inc	1	
p_transmit	1	Fraction of transmitted (only attenuated) neutron rays
DW	1	Global Debye-Waller factor when the 'DW' column is not available. Use 1 if included in F2
nb_atoms	1	Number of sub-unit per unit cell, that is ratio of sigma for chemical formula to sigma per unit cell
d_phi	deg	Angle corresponding to the vertical angular range to focus to, e.g. detector height. 0 for no focusing
p_interact	1	Fraction of events interacting with sample, e.g. 1-p_transmit-p_inc
concentric	only for box, cylinder, sphere	[1] Indicate that this component has a hollow geometry and may contain other components. It should then be duplicated after the inside part
density	g/cm^3	Density of material. rho=density/weight/1e24*N_A.
weight	g/mol	Atomic/molecular weight of material
barns	1	Flag to indicate if F ^2 from 'reflections' is in barns or fm^2 (barns=1 for laz, barns=0 for lau type files).
Strain	ppm	Global relative delta_d_d shift when the 'Strain' column is not available. Use 0 if ideal.
focus_flip	1	Controls the sense of d_phi. If 0 d_phi is measured against the xz-plane. If !=0 d_phi is measured against zy-plane.


Sampling parameters

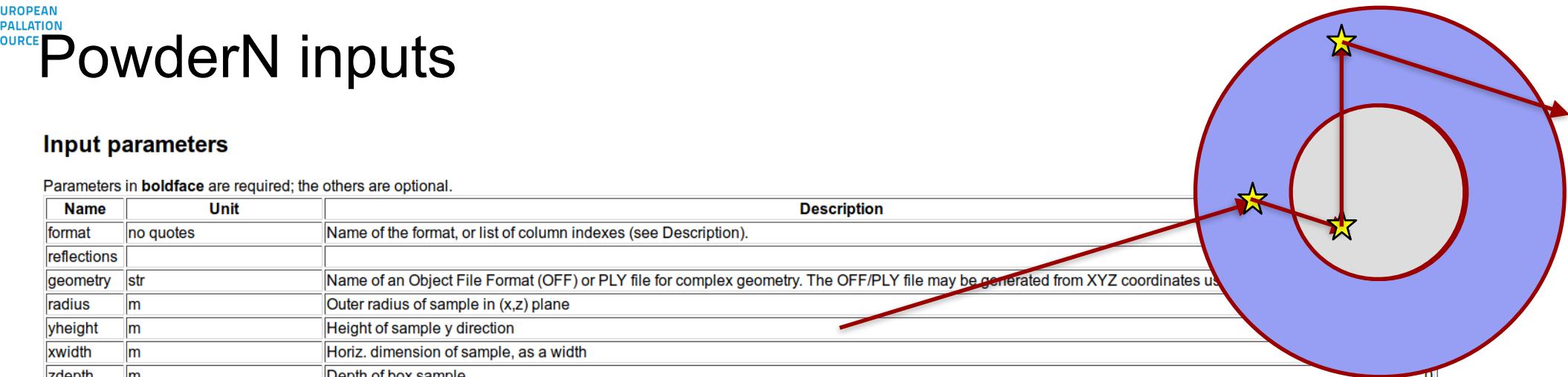


Concentric = Hollow
Powder [cylinder]

Input parametersParameters in **boldface** are re

Name	Unit
format	no quotes
reflections	
geometry	str
radius	m
yheight	m
xwidth	m
zdepth	m
thickness	
pack	1
Vc	AA^3
sigma_abs	barns
sigma_inc	barns
delta_d_d	0/1
p_inc	1
p_transmit	1
DW	1
nb_scans	1
d_phi	deg
p_interact	1
concentric	only for box, cylinder
density	g/cm^3
weight	g/mol
barns	1
Strain	ppm
focus_flip	1

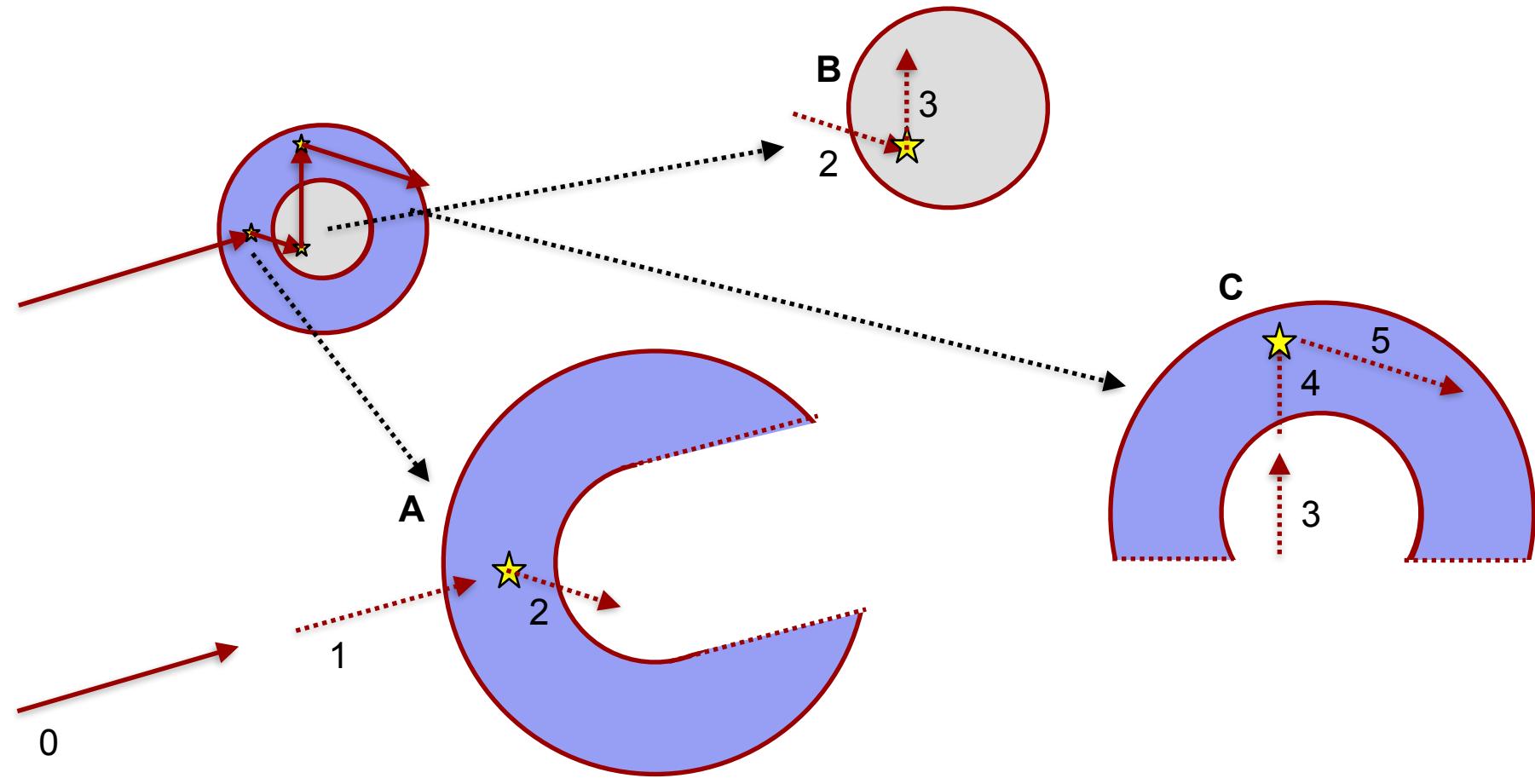**Concentric = Hollow
Powder [cylinder]**


PowderN concentric in code

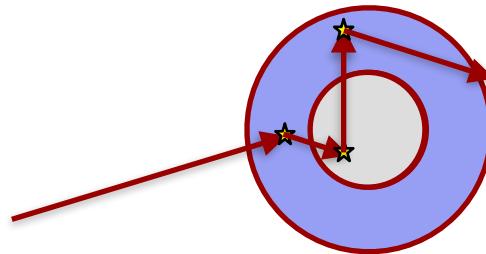



```
COMPONENT Al = PowderN(  
    reflections="Al.laz", radius=0.02, yheight=0.05, thickness=0.005,  
    p_transmit=0.8, concentric=1)  
AT (0, 0, 0) RELATIVE PREVIOUS
```

```
COMPONENT a12co3 = PowderN(  
    reflections=Al2CO3-laz, radius=0.005, yheight=0.05,  
    p_inc=0.15, p_transmit=0.5)  
AT (0, 0, 0) RELATIVE PREVIOUS
```

```
COMPONENT Al_2 = COPY(Al) (concentric=0)  
AT (0,0,0) RELATIVE PREVIOUS
```

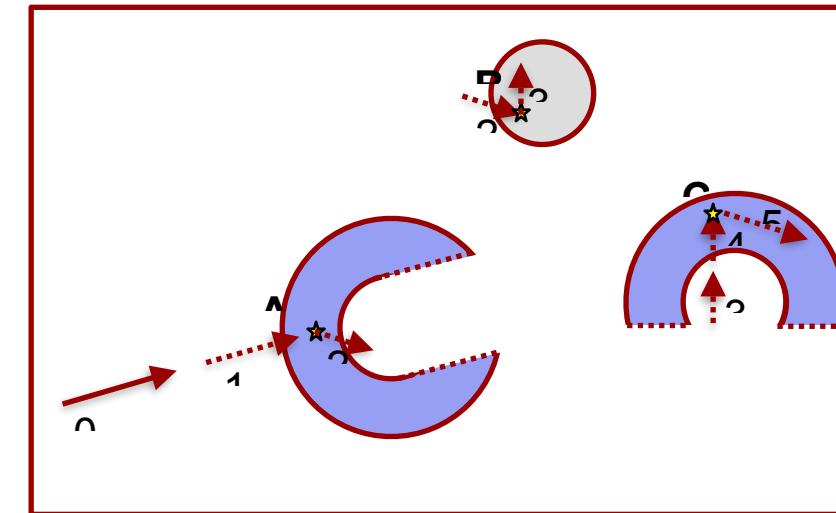

PowderN concentric in code

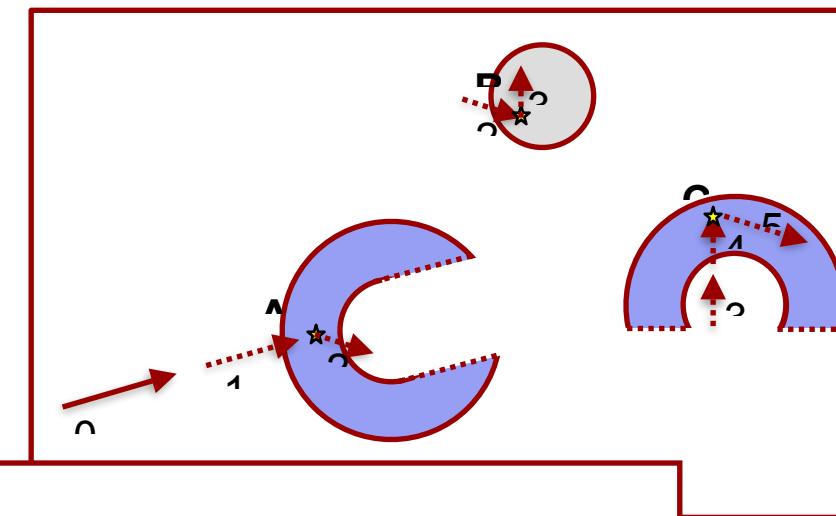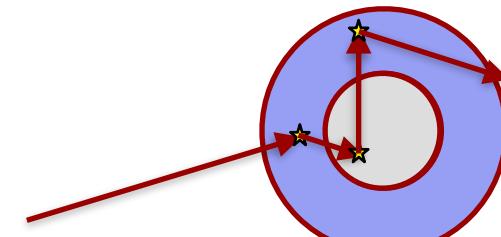

COMPOUND
This is the most basic way to define a material:
AT (0,0,0) * Material (will be) added in “[Supplemental material](#)”

COMPOUND
* A dedicated tutorial is available through the
[McStasScript + Union documentation](#)

AT (0,0,0) * Has most of the material models from “standard” McStas available

COMPOUND
AT (0,0,0) * Is available in McXtrace but not fully ported

AT (0,0,0) [RELATIVE PREVIOUS](#)

PowderN inputs

Input parameters

Parameters in **boldface** are required; the others are optional.

Name	Unit	Description	Default
format	no quotes	Name of the format, or list of column indexes (see Description).	Undefined
reflections			"NULL"
geometry	str	Name of an Object File Format (OFF) or PLY file for complex geometry. The OFF/PLY file may be generated from XYZ coordinates using qhull/powercrust	"NULL"
radius	m	Outer radius of sample in (x,z) plane	0
yheight	m	Height of sample y direction	0
xwidth	m	Horiz. dimension of sample, as a width	0
zdepth	m	Depth of box sample	0
thickness			0
pack	1	Packing factor	1
Vc	AA^3	Volume of unit cell=nb atoms per cell/density of atoms	0
sigma_abs	barns	Absorption cross section per unit cell at 2200 m/s. Use a negative value to unactivate it	0
sigma_inc	barns	Incoherent cross section per unit cell. Use a negative value to unactivate it	0
delta_d_d	0/1	Global relative delta_d_d/d broadening when the 'w' column is not available. Use 0 if ideal.	0
p_inc	1	Fraction of incoherently scattered neutron rays	0.1
p_transmit	1	Fraction of transmitted (only attenuated) neutron rays	0.1
DW	1	Global Debye-Waller factor when the 'DW' column is not available. Use 1 if included in F2	0
nb_atoms	1	Number of sub-unit per unit cell, that is ratio of sigma for chemical formula to sigma per unit cell	1
d_phi	deg	Angle corresponding to the vertical angular range to focus to, e.g. detector height. 0 for no focusing	0
p_interact	1	Fraction of events interacting with sample, e.g. 1-p_transmit-p_inc	0
concentric	only for box, cylinder, sphere	[1] Indicate that this component has a hollow geometry and may contain other components. It should then be duplicated after the inside part	0
density	g/cm^3	Density of material. rho=density/weight/1e24*N_A.	0
weight	g/mol		0
barns	1	s or fm^2 (barns=1 for laz, barns=0 for lau type files).	1
Strain	ppm	column is not available. Use 0 if ideal.	0
focus_flip	1	d against the xz-plane. If !=0 d_phi is measured against zy-plane.	0

Only scatter into a part of
 4π


PowderN inputs

Input parameters

Parameters in **boldface** are required; the others are optional.

Name	Unit	Description	Default
format	no quotes	Name of the format, or list of column indexes (see Description).	Undefined
reflections			"NULL"
geometry	str	Name of an Object File Format (OFF) or PLY file for complex geometry. The OFF/PLY file may be generated from XYZ coordinates using qhull/powercrust	"NULL"
radius	m	Outer radius of sample in (x,z) plane	0
yheight	m	Height of sample y direction	0
xwidth	m	Horiz. dimension of sample, as a width	0
zdepth	m	Depth of box sample	0
thickness			0
pack	1	Packing factor	1
Vc	AA^3	Volume of unit cell=nb atoms per cell/density of atoms	0
sigma_abs	barns	Absorption cross section per unit cell at 2200 m/s. Use a negative value to unactivate it	0
sigma_inc	barns	Incoherent cross section per unit cell. Use a negative value to unactivate it	0
delta_d_d	0/1	Global relative delta_d_d/d broadening when the 'w' column is not available. Use 0 if ideal.	0
p_inc	1	Fraction of incoherently scattered neutron rays	0.1
p_transmit	1	Fraction of transmitted (only attenuated) neutron rays	0.1
DW	1	Global Debye-Waller factor when the 'DW' column is not available. Use 1 if included in F2	0
nb_atoms	1	Number of sub-unit per unit cell, that is ratio of sigma for chemical formula to sigma per unit cell	1
d_phi	deg	Angle corresponding to the vertical angular range to focus to, e.g. detector height. 0 for no focusing	0
p_interact	1	Fraction of events interacting with sample, e.g. 1-p_transmit-p_inc	0
concentric	only for box, cylinder, sphere	[1] Indicate that this component has a hollow geometry and may contain other components. It should then be duplicated after the inside part	0
density	g/cm^3	Density of material. rho=density/weight/1e24*N_A.	0
weight	g/mol		0
barns	1	s or fm^2 (barns=1 for laz, barns=0 for lau type files).	1
Strain	ppm	column is not available. Use 0 if ideal.	0
focus_flip	1	d against the xz-plane. If !=0 d_phi is measured against zy-plane.	0

Only scatter into a part of
 4π

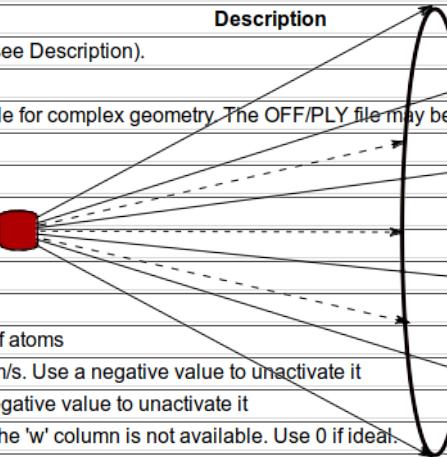

PowderN inputs

Input parameters

Parameters in **boldface** are required; the others are optional.

Name	Unit	Description	Default
format	no quotes	Name of the format, or list of column indexes (see Description).	Undefined
reflections			"NULL"
geometry	str	Name of an Object File Format (OFF) or PLY file for complex geometry. The OFF/PLY file may be generated from XYZ coordinates using qhull/powercrust	"NULL"
radius	m	Outer radius of sample in (x,z) plane	0
yheight	m	Height of sample y direction	0
xwidth	m	Horiz. dimension of sample, as a width	0
zdepth	m	Depth of box sample	0
thickness			0
pack	1	Packing factor	1
Vc	AA^3	Volume of unit cell=nb atoms per cell/density of atoms	0
sigma_abs	barns	Absorption cross section per unit cell at 2200 m/s. Use a negative value to unactivate it	0
sigma_inc	barns	Incoherent cross section per unit cell. Use a negative value to unactivate it	0
delta_d_d	0/1	Global relative delta_d_d/d broadening when the 'w' column is not available. Use 0 if ideal.	0
p_inc	1	Fraction of incoherently scattered neutron rays	0.1
p_transmit	1	Fraction of transmitted (only attenuated) neutron rays	0.1
DW	1	Global Debye-Waller factor when the 'DW' column is not available. Use 1 if included in F2	0
nb_atoms	1	Number of sub-unit per unit cell, that is ratio of sigma for chemical formula to sigma per unit cell	1
d_phi	deg	Angle corresponding to the vertical angular range to focus to, e.g. detector height. 0 for no focusing	0
p_interact	1	Fraction of events interacting with sample, e.g. 1-p_transmit-p_inc	0
concentric	only for box, cylinder, sphere	[1] Indicate that this component has a hollow geometry and may contain other components. It should then be duplicated after the inside part	0
density	g/cm^3	Density of material. rho=density/weight/1e24*N_A.	0
weight	g/mol		0
barns	1	s or fm^2 (barns=1 for laz, barns=0 for lau type files).	1
strain	ppm	column is not available. Use 0 if ideal.	0
focus_flip	1	d against the xz-plane. If !=0 d_phi is measured against zy-plane.	0

Only scatter into a part of
 4π

Reflection files for Single_crystal and PowderN

```
# TITLE *Aluminum-Al-[FM3-M] Miller, H.P.jr.;DuMond, J.W.M.[1942] at 298 K
# CELL 4.049320 4.049320 4.049320 90.000000 90.000000 90.000000
# SPCGRP F M 3 M CUBIC STRUCTURE
# ATOM AL 1 0.000000 0.000000 0.000000
# SCATTERING FACTOR COEFFICIENTS: AL F= 0.345 CM-12
# Reference: Physical Review (1940) 57, 198-206
#
# Physical parameters:
# sigma_coh 1.495 coherent scattering cross section (single atom) in [barn]
# sigma_inc 0.0082 incoherent scattering cross section (single atom)in [barn]
# sigma_abs 0.231 absorption scattering cross section (single atom) in [barn]
# density 2.70 in [g/cm^3]
# weight 26.98 in [g/mol] (single atom)
# multiplicity 4 in [atoms/unit cell]
# Vc 66.4 volume of unit cell in [A^3]
# v_sound 5100 in [m/s]
# v_sound_l 6420 velocity of longitudinal sound in [m/s]
# v_sound_t 3040 velocity of transversal sound in [m/s]
# T_m 933.5 melting temperature in [K]
# T_b 2792.2 boiling temperature in [K]
# At_number 13 atomic number Z
# lattice_a 4.04932 lattice parameter a in [Angs]
#
# Format parameters: Crystallographica format
# column_j 4 multiplicity 'j'
# column_d 5 d-spacing 'd' in [Angs]
# column_F2 7 norm of scattering factor |F|^2 in [fm^2]
# column_h 1
# column_k 2
# column_l 3
#
# h k l Mult. d-space 2Theta F-squared
-1 -1 8 2.338 24.6973 21.3
-1 -1 1 8 2.338 24.6973 21.3
...

```

Lau datafiles

header

+

reflection list

Can be used with
Single_crystal, PowderN,
Isotropic_Sqw

Input for the PowderN

```

# TITLE *Corundum-Al2O3-[R3-CH] Graafsma, H.;Souhassou, M.;Harkem[1998] [corundum saphire:blue, ruby:red]
# CELL 4.757000 4.757000 12.987700 90.000000 90.000000 120.000000
# SPCGRP R -3 C TRIGONAL STRUCTURE
# ATOM AL 1 0.000000 0.000000 0.352110
# ATOM O 1 0.306260 0.306260 0.250000
# SCATTERING FACTOR COEFFICIENTS: AL F= 0.345 CM-12 ; O F= 0.581 CM-12
# Reference: Acta Crystallographica B (1998) 54, 193-195
#
# Physical parameters:
# sigma_coh 15.683 coherent scattering cross section for Al2O3 in [barn]
# sigma_inc 0.0188 incoherent scattering cross section for Al2O3 in [barn]
# sigma_abs 0.4625 absorption scattering cross section for Al2O3 in [barn]
# density 4.05 in [g/cm^3]
# weight 101.96 in [g/mol] for Al2O3
# multiplicity 6 in [Al2O3/unit cell]
# Vc 254.52 volume of unit cell in [A^3]
# T_m 2273 melting temperature in [K]
# T_b 3773 boiling temperature in [K]
# lattice_a 4.757 lattice parameter a in [Angs]
# lattice_c 12.9877 lattice parameter c in [Angs]
# lattice_cc 120 lattice angle gamma in [deg]
#
# Format parameters: Lazy format <http://icsd.ill.fr>
# column_j 17 multiplicity 'j'
# column_d 6 d-spacing 'd' in [Angs]
# column_F 13 norm of scattering factor |F| in [barn]
# column_h 1
# column_k 2
# column_l 3
#
# H K L THETA 2THETA D VALUE 1/D**2 SIN2*1000 H K L INTENSITY /F(HKL)/ A(HKL) B(HKL) PHA.ANG. MULT LPG
1 0 1 6.35 12.71 4.5175 0.0490 12.25 1 0 1 367.0 4.1 -4.08 0.00 180.00 6 82.14
0 0 3 7.10 14.20 4.0467 0.0611 15.27 0 0 3 110.0 4.3 4.32 0.00 0.00 2 66.01
0 1 2 7.57 15.13 3.7972 0.0694 17.34 0 1 2 10.9 0.8 0.84 0.00 0.00 6 58.18
...

```

Laz + Lau datafiles

header

+

reflection list

Can be used with PowderN,
Isotropic_Sqw

What's already there?

Single Crystals

Al.lau
 Al203_sapphire.lau
 B4C.lau
 BCC_fake.lau
 BeO.lau
 C60.lau
 C_diamond.lau
 C_graphite.lau
 CaF2.lau
 Ge.lau
 In.lau
 PPase_D_P1.lau
 Rubredoxin.lau
 SiO2_cristobalite.lau
 SiO2_quartz.a.lau
 SiO2_quartzb.lau
 TiO2_rutile.lau
 YBaCuO.lau
 adrenaline.lau
 aspirin.lau
 leucine.lau

Powders


Ag.laz
 Al.laz
 Al203_sapphire.laz
 Au.laz
 B4C.laz
 Ba.laz
 Be.laz
 BeO.laz
 Bi.laz
 C_diamond.laz
 C_graphite.laz
 Cd.laz
 Cr.laz
 Cs.laz
 Cu.laz
 Cu2MnAl.laz
 Fe.laz
 Fe_Gamma.laz
 Ga.laz
 Gd.laz
 Ge.laz
 H2O_ice_1h.laz
 He4_hcp.laz
 Hg.laz
 I2.laz
 K.laz
 Li.laz
 LiF.laz
 Mo.laz
 Na2Ca3Al2F14.laz
 Nb.laz
 Ni.laz
 Pb.laz
 Pd.laz
 Pt.laz
 Rb.laz
 Se_alpha.laz
 Se_beta.laz
 Si.laz
 SiO2_quartz.a.laz
 SiO2_quartzb.laz
 Sn_alpha.laz
 Sn_beta.laz
 Ti.laz
 Tl.laz
 UO2.laz
 V.laz
 Y2O3.laz
 Y3Fe5O12_YIG.laz
 Zn.laz
 Zr.laz
 duplex.laz

Where to get these files...

- | \$MCSTAS/data - use mcdoc
- | - Or make your own via
 - | Finding a CIF file for the given structure
 - | e.g. from
 - | ICSD <http://icsd.fiz-karlsruhe.de> or
 - | COD <<http://crystallography.net>>
 - | In the current McStas, PowderN, Single_crystal “understands” CIF input - applies
 - | cif2hkl behind the scenes
 - | cif2hkl is available in your McStas installation

More detailed info available here:

https://github.com/McStasMcXtrace/McCode/wiki/Generating-input-files-for-Single_crystal-and-PowderN

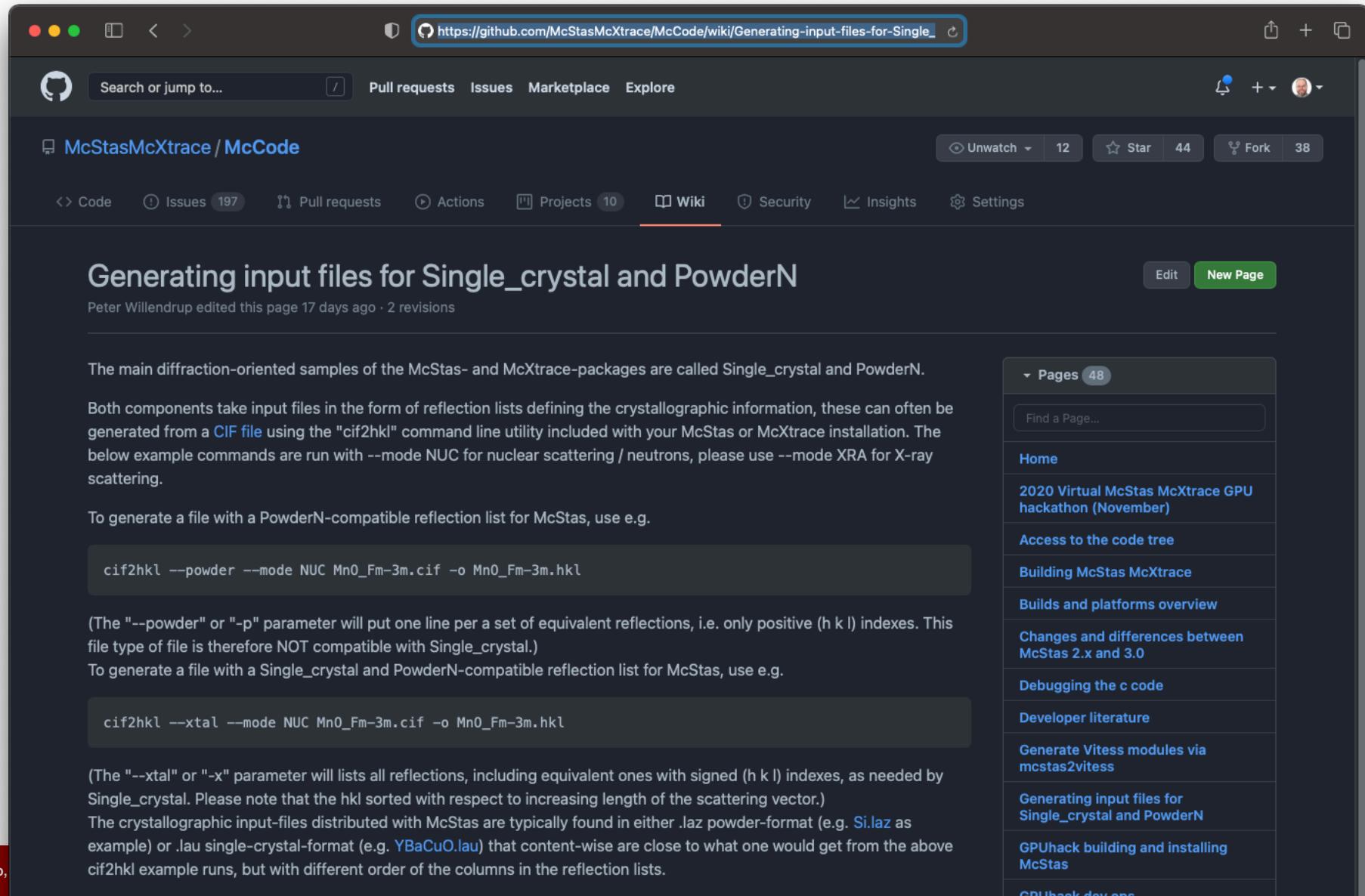

The screenshot shows a GitHub repository page for 'McStasMcXtrace / McCode'. The 'Wiki' tab is selected. The main content is a page titled 'Generating input files for Single_crystal and PowderN'. It was last edited by Peter Willendrup 17 days ago. The page contains text and code snippets explaining how to generate reflection lists for McStas and McXtrace.

Generating input files for Single_crystal and PowderN

Peter Willendrup edited this page 17 days ago · 2 revisions

The main diffraction-oriented samples of the McStas- and McXtrace-packages are called Single_crystal and PowderN.

Both components take input files in the form of reflection lists defining the crystallographic information, these can often be generated from a CIF file using the "cif2hkl" command line utility included with your McStas or McXtrace installation. The below example commands are run with --mode NUC for nuclear scattering / neutrons, please use --mode XRA for X-ray scattering.

To generate a file with a PowderN-compatible reflection list for McStas, use e.g.

```
cif2hkl --powder --mode NUC Mn0_Fm-3m.cif -o Mn0_Fm-3m.hkl
```

(The "--powder" or "-p" parameter will put one line per a set of equivalent reflections, i.e. only positive (h k l) indexes. This file type of file is therefore NOT compatible with Single_crystal.)

To generate a file with a Single_crystal and PowderN-compatible reflection list for McStas, use e.g.

```
cif2hkl --xtal --mode NUC Mn0_Fm-3m.cif -o Mn0_Fm-3m.hkl
```

(The "--xtal" or "-x" parameter will lists all reflections, including equivalent ones with signed (h k l) indexes, as needed by Single_crystal. Please note that the hkl sorted with respect to increasing length of the scattering vector.)

The crystallographic input-files distributed with McStas are typically found in either .laz powder-format (e.g. Si.laz as example) or .lau single-crystal-format (e.g. YBaCuO.lau) that content-wise are close to what one would get from the above cif2hkl example runs, but with different order of the columns in the reflection lists.

Pages 48

- Find a Page...
- Home
- 2020 Virtual McStas McXtrace GPU hackathon (November)
- Access to the code tree
- Building McStas McXtrace
- Builds and platforms overview
- Changes and differences between McStas 2.x and 3.0
- Debugging the c code
- Developer literature
- Generate Vitess modules via mcstas2vitess
- Generating input files for Single_crystal and PowderN
- GPUhack building and installing McStas
- GPUhack dev doc

McStas Data formats

McStas uses a number of input data formats, all text based.

- **Lau** (Single_crystal, PowderN, Isotropic_Sqw)
- **Laz** (PowderN, Isotropic_Sqw)
- **Sqw** (Isotropic_Sqw)
- **qSq** (Isotropic_Sqw) – mostly for liquids
- **PDB** (SANS_PDB)
- **ncmat** (NCrystal)

OFF/PLY any geometry

An additional complex geometry enables to use any point set to describe the material volume (geomview OFF file).

