



Erik Knudsen

# McStas samples for diffraction

# Agenda

- General concepts (reminder)
- Focus on diffraction-oriented samples:
  - Incoherent scatterer
  - Powder
  - Single Crystal
  - File formats

# Incoherent Sample

## [Incoherent.comp](#)

Parameters in **boldface** are required; the others are optional.

An incoherent scatterer with various sample shape options

Name	Unit	Description	Default
geometry	str	Name of an Object File Format (OFF) or PLY file for complex geometry. The OFF/PLY file may be generated from XYZ coordinates using qhull/powercrust	0
radius	m	Outer radius of sample in (x,z) plane	0
xwidth	m	Horiz. dimension of sample (bounding box if off file), as a width	0
yheight	m	Vert. dimension of sample (bounding box if off file), as a height. A sphere shape is used when 0 and radius is set	0
zdepth	m	Depth of sample (bounding box if off file)	0
thickness	m	Thickness of hollow sample	0
target_x	-		0
target_y	m	position of target to focus at	0
target_z	-		0
focus_r	m	Radius of disk containing target. Use 0 for full space	0
focus_xw	m	horiz. dimension of a rectangular area	0
focus_yh	m	vert. dimension of a rectangular area	0
focus_aw	deg	horiz. angular dimension of a rectangular area	0
focus_ah	deg	vert. angular dimension of a rectangular area	0
target_index	1	Relative index of component to focus at, e.g. next is +1	0
pack	1	Packing factor	1
p_interact	1	MC Probability for scattering the ray; otherwise transmit	1
f_QE	1	Fraction of quasielastic scattering (rest is elastic)	0
gamma	1	Lorentzian width of quasielastic broadening (HWHM)	0
sigma_abs	barns	Absorption cross section pr. unit cell at 2200 m/s	5.08
sigma_inc	barns	Incoherent scattering cross section pr. unit cell	5.08
Vc	AA <sup>3</sup>	Unit cell volume	13.827
concentric	1	Indicate that this component has a hollow geometry and may contain other components. It should then be duplicated after the inside part (only for box, cylinder, sphere)	0
order	-	Limit multiple scattering up to given order	0

# Incoherent Sample

## Incoherent.comp

Parameters in **boldface** are required; the others are optional.

Name	Unit	Description	Default
geometry	str	Name of an Object File Format (OFF) or PLY file for complex geometry. The OFF/PLY file may be generated from XYZ coordinates using qhull/powercrust	0
radius	m	Outer radius of sample in (x,z) plane	0
xwidth	m	Horiz. dimension of sample (bounding box if off file), as a width	0
yheight	m	Vert. dimension of sample (bounding box if off file), as a height. A sphere shape is used when 0 and radius is set	0
zdepth	m	Depth of sample (bounding box if off file)	0
thickness	m	Thickness of hollow sample	0
target_x	-		0
target_y	m	position of target to focus at	0
target_z	-		0
focus_r	m	Radius of disk containing target. Use 0 for full space	0
focus_xw	m	horiz. dimension of a rectangular area	0
focus_yh	m	vert. dimension of a rectangular area	
focus_aw	deg	horiz. angular dimension of a rectangular area	
focus_ah	deg	vert. angular dimension of a rectangular area	
target_index	1	Relative index of component in stack	
pack	1	Packing factor	
p_interact	1	MC Probability for scattering	
I_QE	1	Fraction of quasielastic intensity	
gamma	1	Lorentzian width of quasielastic peak	
sigma_abs	barns	Absorption cross section	
sigma_inc	barns	Incoherent scattering cross section	
Vc	AA <sup>3</sup>	Unit cell volume	13.827
concentric	1	Indicate that this component has a hollow geometry and may contain other components. It should then be duplicated after the inside part (only for box, cylinder, sphere)	0
order	-	Limit multiple scattering up to given order	0

An incoherent scatterer with various sample shape options


**p\_interact** = The probability that a given ray will interact (scatter) with a sample. Does not directly affect the intensity – only the quality of the resulting Monte Carlo estimate.

# Crystalline Samples





# From single crystal / crystallites to powder....


· Single crystal

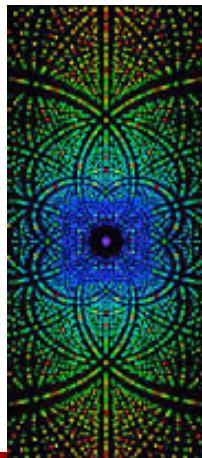


· Polycrystal with a little disorder,  
i.e. a *preferred orientation, texture*

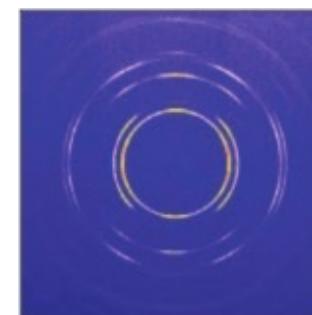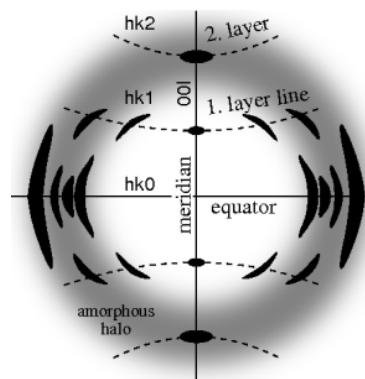

· Powder with complete disorder

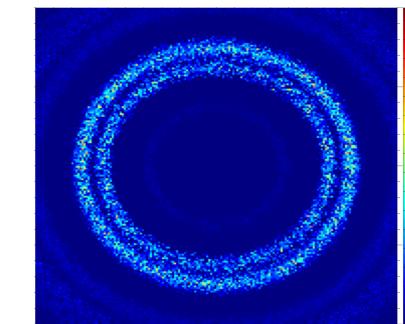

# Single Crystal model

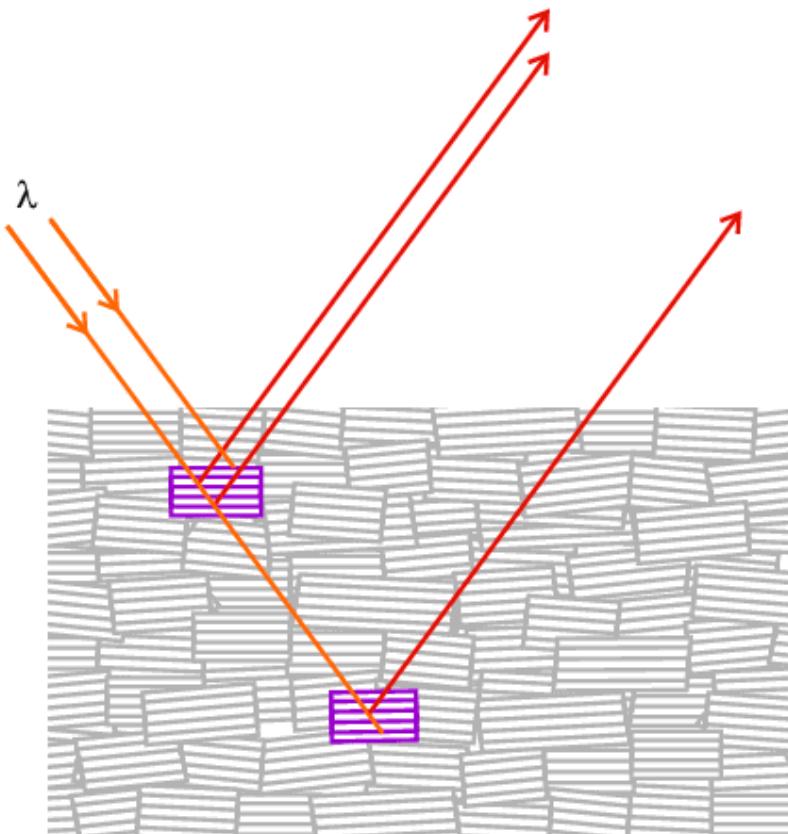

- *Models ideally imperfect crystal:*
- *Peak broadening is dominated by mosaic*  
=> not appropriate for highly perfect crystals.
- “Space group symmetries” considered out of scope.
  - Optimized for small unit cells but...
  - No assumption about unit cell, orientation etc.
  - Gaussian model for mosaic (see later)

Courtesy University College London

# Single Crystal

## Input parameters

Parameters in **boldface** are required; the others are optional.

Name	Unit	Description	Default
mosaic_AB	arc_minutes, arc_minutes,1, 1, 1, 1, 1, 1	In Plane mosaic rotation and plane vectors (anisotropic), mosaic_A, mosaic_B, A_h,A_k,A_l, B_h,B_k,B_l. Puts the crystal in the in-plane mosaic state. Vectors A and B define plane in which the crystal roation is defined, and mosaic_A, mosaic_B, denotes the resp. mosaicities (gaussian RMS) with respect to the two reflections chosen by A and B (Miller indices).	Mosaic_AB_Undefined
reflections	string	File name containing structure factors of reflections. Use empty ("") or NULL for incoherent scattering only	0
geometry	str	Name of an Object File Format (OFF) or PLY file for complex geometry. The OFF/PLY file may be generated from XYZ coordinates using qhull/powercrust	0
xwidth	m	Width of crystal	0
yheight	m	Height of crystal	0
zdepth	no extinction simulated	[m] Depth of crystal	0
radius	m	Outer radius of sample in (x,z) plane	0
delta_d_d	1	Lattice spacing variance, gaussian RMS	1e-4
mosaic	arc minutes	Crystal mosaic (isotropic), gaussian RMS. Puts the crystal in the isotropic mosaic model state, thus disregarding other mosaicity parameters.	-1
mosaic_a	arc minutes	Horizontal (rotation around lattice vector a) mosaic (anisotropic), gaussian RMS. Put the crystal in the anisotropic crystal vector state. I.e. model mosaicity through rotation around the crystal lattice vectors. Has precedence over in-plane mosaic model.	-1
mosaic_b	arc minutes	Vertical (rotation around lattice vector b) mosaic (anisotropic), gaussian RMS.	-1
mosaic_c	arc minutes	Out-of-plane (Rotation around lattice vector c) mosaic (anisotropic), gaussian RMS	-1
recip_cell	1	Choice of direct/reciprocal (0/1) unit cell definition	0
barns	1	Flag to indicate if $ F ^2$ from 'reflections' is in barns or fm <sup>2</sup> . barns=1 for laz and isotropic constant elastic scattering (reflections=NULL), barns=0 for lau type files	0
ax	AA or AA <sup>-1</sup>	Coordinates of first (direct/recip) unit cell vector	0
ay			0
az			0
bx	AA or AA <sup>-1</sup>	Coordinates of second (direct/recip) unit cell vector	0
by			0
bz			0
cx	AA or AA <sup>-1</sup>	Coordinates of third (direct/recip) unit cell vector	0
cy			0
cz			0
p_transmit	1	Monte Carlo probability for neutrons to be transmitted without any scattering. Used to improve statistics from weak reflections	-1
sigma_abs	barns	Absorption cross-section per unit cell at 2200 m/s	0
sigma_inc	barns	Incoherent scattering cross-section per unit cell Use -1 to unactivate	0
aa	deg	Unit cell angles alpha, beta and gamma. Then uses norms of vectors a,b and c as lattice parameters	0
bb	deg	Beta angle	0
cc	deg	Gamma angle	0
order	0: all, 1: first, 2: second, ...	[1] Limit multiple scattering up to given order	0
RX	m	Radius of horizontal along X lattice curvature. flat for 0	0
RY	m	Radius of vertical lattice curvature. flat for 0	0
RZ	m	Radius of horizontal along Z lattice curvature. flat for 0	0
powder	1	Flag to indicate powder mode, for simulation of Debye-Scherrer cones via random crystallite orientation. A powder texture can be approximated with 0	0
PG	1	Flag to indicate "Pyrolytic Graphite" mode, only meaningful with choice of Graphite.lau, models PG crystal. A powder texture can be approximated with 0	0

# Single Crystal

## Input parameters

Parameters in **boldface** are required; the others are optional.

Name	Unit	Description	Default
mosaic_AB	arc_minutes, arc_minutes, 1, 1, 1, 1, 1	In Plane mosaic rotation and plane vectors (anisotropic), mosaic_A, mosaic_B, A_h,A_k,A_l, B_h,B_k,B_l. Puts the crystal in the in-plane mosaic state. Vectors A and B define plane in which the crystal roation is defined, and mosaic_A, mosaic_B, denotes the resp. mosaicities (gaussian RMS) with respect to the two reflections chosen by A and B (Miller indices).	Mosaic_AB_Undefined
reflections	string	File name containing structure factors of reflections. Use empty ("") or NULL for incoherent scattering only	0
geometry	str	Name of an Object File Format (OFF) or PLY file for complex geometry. The OFF/PLY file may be generated from XYZ coordinates using qhull/powercrust	0
xwidth	m	Width of crystal	0
yheight	m	Height of crystal	0
zdepth	no extinction simulated	[m] Depth of crystal	0
radius	m	Outer radius of sample in (x,z) plane	0
delta_d_d	1	Lattice spacing variance, gaussian RMS	1e-4
mosaic	arc minutes	Crystal mosaic (isotropic), gaussian RMS. Puts the crystal in the isotropic mosaic model state, thus disregarding other mosaicity parameters.	-1
mosaic_a	arc minutes	Horizontal (rotation around lattice vector a) mosaic (anisotropic), gaussian RMS. Put the crystal in the anisotropic crystal vector state. I.e. model mosaicity through rotation around the crystal lattice vectors. Has precedence over in-plane mosaic model.	-1
mosaic_b	arc minutes	Vertical (rotation around lattice vector b) mosaic (anisotropic), gaussian RMS.	-1
mosaic_c	arc minutes	Out-of-plane (Rotation around lattice vector c) mosaic (anisotropic), gaussian RMS	-1
recip_cell	1	Choice of direct/reciprocal (0/1) unit cell definition	0
barns	1	Flag to indicate if $ F ^2$ from 'reflections' is in barns or fm^2. barns=1 for laz and isotropic constant elastic scattering (reflections=NULL), barns=0 for lau type files	0
ax	AA or AA^-1	Coordinates of first (direct/recip) unit cell vector	0
ay			0
az			0
bx	AA or AA^-1	Coordinates of second (direct/recip) unit cell vector	0
by			0
bz			0
cx	AA or AA^-1		0
cy			0
cz			0
p_transm	1		-1
sigma_pos	barns		0
sigm_a_inc	barns	Incoherent scattering cross-section per unit cell Use -1 to deactivate	0
aa	deg	Unit cell angles alpha, beta and gamma. Then uses norms of vectors a,b and c as lattice parameters	0
bb	deg	Beta angle	0
cc	deg	Gamma angle	0
order	0: all, 1: first, 2: second, ...	[1] Limit multiple scattering up to given order	0
RX	m	Radius of horizontal along X lattice curvature. flat for 0	0
RY	m	Radius of vertical lattice curvature. flat for 0	0
RZ	m	Radius of horizontal along Z lattice curvature. flat for 0	0
powder	1	Flag to indicate powder mode, for simulation of Debye-Scherrer cones via random crystallite orientation. A powder texture can be approximated with 0	0
PG	1	Flag to indicate "Pyrolytic Graphite" mode, only meaningful with choice of Graphite.lau, models PG crystal. A powder texture can be approximated with 0	0

Unit cell vectors (ax,ay,az, bx, by,bz, cx,cy, cz). Encodes the crystal orientation relative to the Component's orientation

# Single Crystal


## Input parameters

Parameters in **boldface** are required; the others are optional.

Name	Unit	Description	Default
<b>mosaic_AB</b>	arc_minutes, arc_minutes,1, 1, 1	In Plane mosaic rotation and plane vectors (anisotropic), mosaic_A, mosaic_B, A_h.A_k,A_l, B_h.B_k,B_l. Puts the crystal in the in-plane mosaic state. Vectors A and B define plane in which the crystal roation is defined, and mosaic_A, mosaic_B, denotes the resp. mosaicities (gaussian RMS) with respect to the two reflections chosen by A and B (Miller indices).	Mosaic_AB_Undefined
<b>reflections</b>	string	File name containing structure factors of reflections. Use empty ("") or NULL for incoherent scattering only	0
<b>g_mosaic</b>	str	Name of an Object File Format (OFF) or PLY file for complex geometry. The OFF/PLY file may be generated from XYZ coordinates using qhull/powercrust	0
xwidth	m	Width of crystal	0
yheight	m	Height of crystal	0
zdepth	no extinction simulated	[m] Depth of crystal	0
radius	m	Outer radius of sample in (x,z) plane	0
<b>lens_d</b>	1		1e-4
<b>mosaic</b>	arc minutes		-1
<b>mosaic_a</b>	arc minutes		-1
<b>mosaic_b</b>	arc minutes		-1
<b>mosaic_c</b>	arc minutes		-1
<b>recip_cell</b>	1		0
<b>barns</b>	1	Flag to indicate if $ F ^2$ from 'reflections' is in barns or fm $^2$ . barns=1 for laz and isotropic constant elastic scattering (reflections=NULL), barns=0 for lau type files	0
<b>cx</b>	AA or AA $^{-1}$	Coordinates of first (direct/recip) unit cell vector	0
<b>ay</b>			0
<b>az</b>			0
<b>bx</b>	AA or AA $^{-1}$	Coordinates of second (direct/recip) unit cell vector	0
<b>by</b>			0
<b>bz</b>			0
<b>cx</b>	AA or AA $^{-1}$	Coordinates of third (direct/recip) unit cell vector	0
<b>cy</b>			0
<b>cz</b>			0
<b>p_transmit</b>	1	Monte Carlo probability for neutrons to be transmitted without any scattering. Used to improve statistics from weak reflections	-1
<b>sigma_abs</b>	barns	Absorption cross-section per unit cell at 2200 m/s	0
<b>sigma_inc</b>	barns	Incoherent scattering cross-section per unit cell Use -1 to unactivate	0
<b>aa</b>	deg	Unit cell angles alpha, beta and gamma. Then uses norms of vectors a,b and c as lattice parameters	0
<b>bb</b>	deg	Beta angle	0
<b>cc</b>	deg	Gamma angle	0
<b>order</b>	0, all, 1: first, 2: second, ...	[1] Limit multiple scattering up to given order	0
<b>RX</b>	m	Radius of horizontal along X lattice curvature. flat for 0	0
<b>RY</b>	m	Radius of vertical lattice curvature. flat for 0	0
<b>RZ</b>	m	Radius of horizontal along Z lattice curvature. flat for 0	0
<b>powder</b>	1	Flag to indicate powder mode, for simulation of Debye-Scherrer cones via random crystallite orientation. A powder texture can be approximated with 0	0
<b>PG</b>	1	Flag to indicate "Pyrolytic Graphite" mode, only meaningful with choice of Graphite.lau, models PG crystal. A powder texture can be approximated with 0	0

Mosaicity parameters: mosaic, mosaic\_a, mosaic\_b, mosaic\_c, mosaic\_AB


# Single Crystal mosaicity

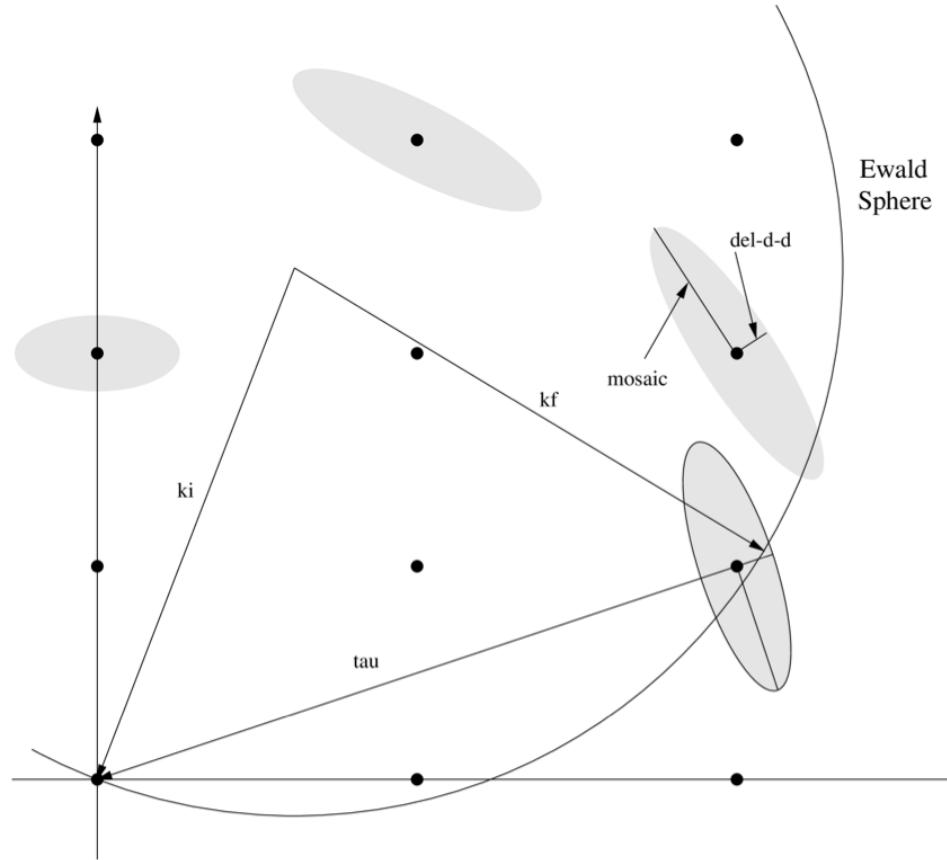

# Single Crystal mosaicity

Three models:



- Isotropic Mosaicity = Gaussian distribution of rotation angle, with given standard deviation, around a random vector.
- Anistropic mosaicity = Gaussian distribution of angles around each of the unit cell axis.

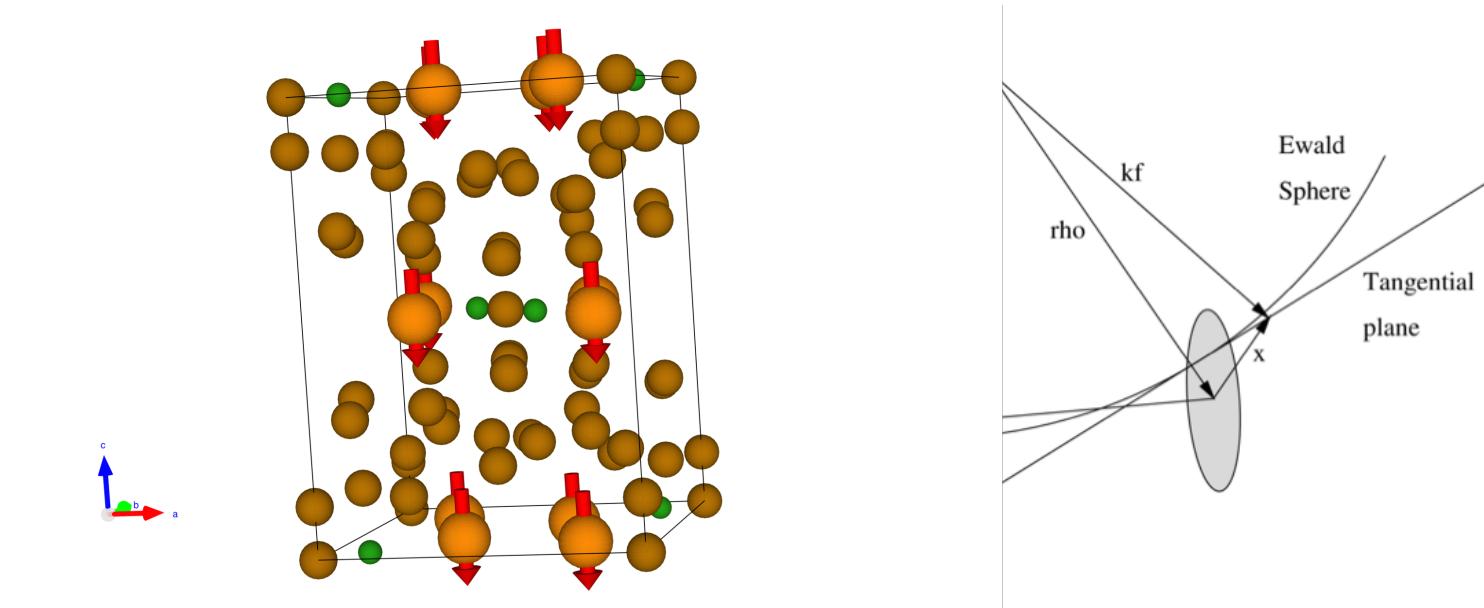

# Single crystal mosaicity

Three models:

- Isotropic Mosaicity = Gaussian distribution of rotation angle, with given standard deviation, around a random vector.
- Anistropic mosaicity = Gaussian distribution of angles around each of the unit cell axis.
- “Rocking curve mosaic” = Gaussian distr. around axis determined by width of two given Bragg peaks. Specify  $\sigma$ , and  $hkl$  for two peaks.

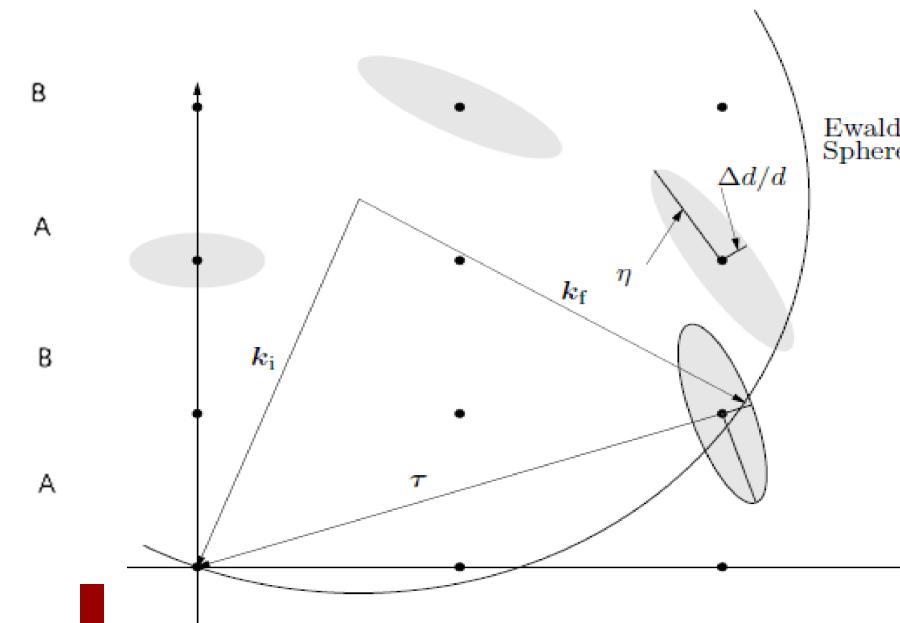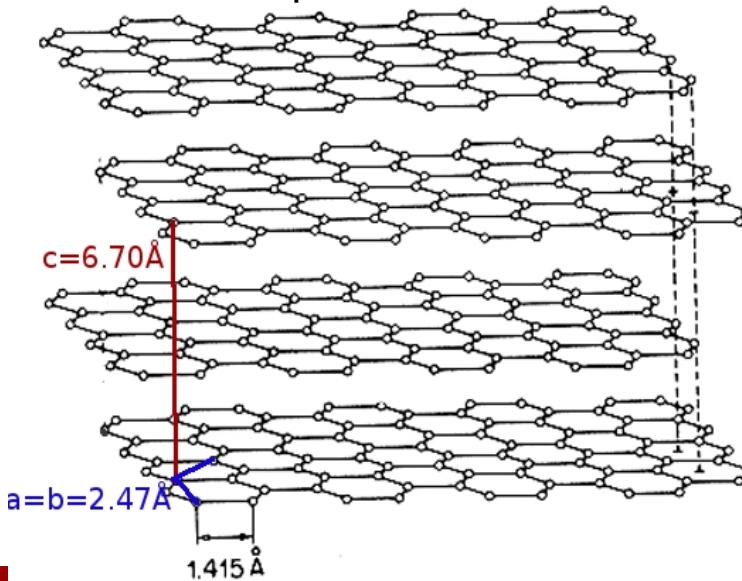

## 9.2 PowderN

Among the most common materials measured with neutron scattering techniques are powders, which are composed of a large number of tiny single crystals. The scattering intensity, by means of constructive interference and averaging over all crystallites, appears as well defined rings around the out-going direct beam direction. These rings are also present in liquids, but are much smoother, as we shall see in next section.

Let us consider a material of density  $\rho$  made of  $N$  scattering units each with a unit cell volume  $V_0$  arranged with lattice spacings  $d_Q=2\pi/Q$  associated to structure factors  $F(Q)$  with equivalent reflection multiplicities  $j_Q$ . These structure factors characterize the efficiency of the reflection with momentum exchange  $Q$ . Following Squires [14] the scattering probability for an incoming neutron with wavelength  $\lambda$  penetrating along a distance  $x$  into the material is about  $1-\exp(-\rho\sigma_{cone}x)$  where  $\sigma_{cone}$  is the so-called coherent elastic cross-section of the ring,

$$\sigma_{cone} = \frac{N \pi \lambda^2 j_Q |F(Q)|^2}{V_0 Q} .$$

This relation is only valid under certain conditions, among which  $d_Q > \lambda/2$ , from the Bragg law. As the possible lattice spacings  $d$  in the material can not exceed a maximum value (for instance the inter-atomic distance), it appears that as the neutron wavelength increases, the number of visible rings in the diffractogram will decrease, until no more scattering is possible above the so-called *Bragg edge*, where materials become transparent to neutrons (except for absorption and incoherent scattering). This is why most of the diffractometers use thermal and hot neutrons. Cold neutrons can only scatter on large distance arrangements in materials, *e.g.* in larger molecules and proteins.

# PowderN inputs

## Input parameters

Parameters in **boldface** are required; the others are optional.

Name	Unit	Description	Default
format	no quotes	Name of the format, or list of column indexes (see Description).	Undefined
reflections			"NULL"
geometry	str	Name of an Object File Format (OFF) or PLY file for complex geometry. The OFF/PLY file may be generated from XYZ coordinates using qhull/powercrust	"NULL"
radius	m	Outer radius of sample in (x,z) plane	0
yheight	m	Height of sample y direction	0
xwidth	m	Horiz. dimension of sample, as a width	0
zdepth	m	Depth of box sample	0
thickness			0
pack	1	Packing factor	1
Vc	AA^3	Volume of unit cell=nb atoms per cell/density of atoms	0
sigma_abs	barns	Absorption cross section per unit cell at 2200 m/s. Use a negative value to unactivate it	0
sigma_inc	barns	Incoherent cross section per unit cell. Use a negative value to unactivate it	0
delta_d_d	0/1	Global relative delta_d_d/d broadening when the 'w' column is not available. Use 0 if ideal.	0
p_inc	1	Fraction of incoherently scattered neutron rays	0.1
p_transmit	1	Fraction of transmitted (only attenuated) neutron rays	0.1
DW	1	Global Debye-Waller factor when the 'DW' column is not available. Use 1 if included in F2	0
nb_atoms	1	Number of sub-unit per unit cell, that is ratio of sigma for chemical formula to sigma per unit cell	1
d_phi	deg	Angle corresponding to the vertical angular range to focus to, e.g. detector height. 0 for no focusing	0
p_interact	1	Fraction of events interacting with sample, e.g. 1-p_transmit-p_inc	0
concentric	only for box, cylinder, sphere	[1] Indicate that this component has a hollow geometry and may contain other components. It should then be duplicated after the inside part	0
density	g/cm^3	Density of material. rho=density/weight/1e24*N_A.	0
weight	g/mol	Atomic/molecular weight of material	0
barns	1	Flag to indicate if  F ^2 from 'reflections' is in barns or fm^2 (barns=1 for laz, barns=0 for lau type files).	1
Strain	ppm	Global relative delta_d_d/d shift when the 'Strain' column is not available. Use 0 if ideal.	0
focus_flip	1	Controls the sense of d_phi. If 0 d_phi is measured against the xz-plane. If !=0 d_phi is measured against zy-plane.	0



# PowderN inputs

## Input parameters

Parameters in **boldface** are required; the others are optional.

Name	Unit	Description	Default
format	no quotes	Name of the format, or list of column indexes (see Description).	Undefined
reflections			"NULL"
geometry	str	Name of an Object File Format (OFF) or PLY file for complex geometry. The OFF/PLY file may be generated from XYZ coordinates using qhull/powercrust	"NULL"
radius	m	Outer radius of sample in (x,z) plane	0
yheight	m	Height of sample y direction	0
xwidth	m	Horiz. dimension of sample, as a width	0
zdepth	m	Depth of box sample	0
thickness			0
pack	1	Packing factor	1
Vc	AA^3	Volume of unit cell=nb atoms per cell/density of atoms	0
sigma_abs	barns	Absorption cross section per unit cell at 2200 m/s. Use a negative value to unactivate it	0
sigma_inc	barns	Incoherent cross section per unit cell. Use a negative value to unactivate it	0
delta_d_d	0/1		0
p_inc	1		0.1
p_transmit	1	Fraction of transmitted (only attenuated) neutron rays	0.1
DW	1	Global Debye-Waller factor when the 'DW' column is not available. Use 1 if included in F2	0
nb_atoms	1	Number of sub-unit per unit cell, that is ratio of sigma for chemical formula to sigma per unit cell	1
d_phi	deg	Angle corresponding to the vertical angular range to focus to, e.g. detector height. 0 for no focusing	0
p_interact	1	Fraction of events interacting with sample, e.g. 1-p_transmit-p_inc	0
concentric	only for box, cylinder, sphere	[1] Indicate that this component has a hollow geometry and may contain other components. It should then be duplicated after the inside part	0
density	g/cm^3	Density of material. rho=density/weight/1e24*N_A.	0
weight	g/mol	Atomic/molecular weight of material	0
barns	1	Flag to indicate if  F ^2 from 'reflections' is in barns or fm^2 (barns=1 for laz, barns=0 for lau type files).	1
Strain	ppm	Global relative delta_d_d shift when the 'Strain' column is not available. Use 0 if ideal.	0
focus_flip	1	Controls the sense of d_phi. If 0 d_phi is measured against the xz-plane. If !=0 d_phi is measured against zy-plane.	0

### Sampling parameters

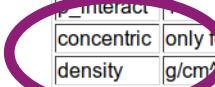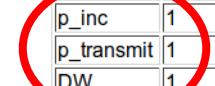

Concentric = Hollow  
Powder [cylinder]


# PowderN inputs

## Input parameters

Parameters in **boldface** are required; the others are optional.


Name	Unit	Description
format	no quotes	Name of the format, or list of column indexes (see Description).
reflections		
geometry	str	Name of an Object File Format (OFF) or PLY file for complex geometry. The OFF/PLY file may be generated from XYZ coordinates us
radius	m	Outer radius of sample in (x,z) plane
yheight	m	Height of sample y direction
xwidth	m	Horiz. dimension of sample, as a width
zdepth	m	Depth of box sample
thickness		
pack	1	Packing factor
Vc	AA^3	Volume of unit cell=nb atoms per cell/density of atoms
sigma_abs	barns	Absorption cross section per unit cell at 2200 m/s. Use a negative value to unactivate it
sigma_inc	barns	Incoherent cross section per unit cell. Use a negative value to unactivate it
delta_d_d	0/1	In the 'w' column is not available. Use 0 if ideal.
p_inc	1	
p_transmit	1	Fraction of transmitted (only attenuated) neutron rays
DW	1	Global Debye-Waller factor when the 'DW' column is not available. Use 1 if included in F2
nb_atoms	1	Number of sub-unit per unit cell, that is ratio of sigma for chemical formula to sigma per unit cell
d_phi	deg	Angle corresponding to the vertical angular range to focus to, e.g. detector height. 0 for no focusing
p_interact	1	Fraction of events interacting with sample, e.g. 1-p_transmit-p_inc
concentric	only for box, cylinder, sphere	[1] Indicate that this component has a hollow geometry and may contain other components. It should then be duplicated after the inside part
density	g/cm^3	Density of material. rho=density/weight/1e24*N_A.
weight	g/mol	Atomic/molecular weight of material
barns	1	Flag to indicate if  F ^2 from 'reflections' is in barns or fm^2 (barns=1 for laz, barns=0 for lau type files).
Strain	ppm	Global relative delta_d_d shift when the 'Strain' column is not available. Use 0 if ideal.
focus_flip	1	Controls the sense of d_phi. If 0 d_phi is measured against the xz-plane. If !=0 d_phi is measured against zy-plane.

Concentric = Hollow  
Powder [cylinder]

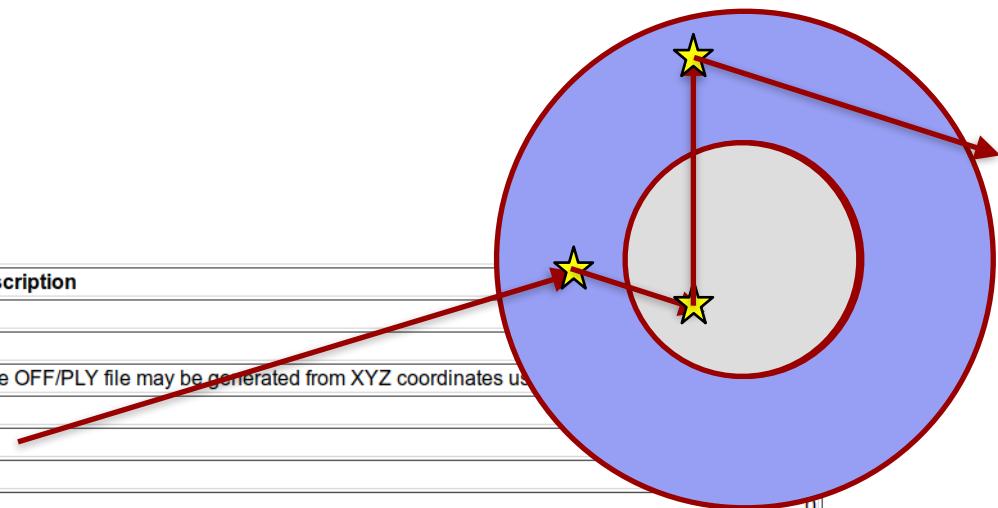


**Input parameters**Parameters in **boldface** are re

Name	Unit
format	no quotes
reflections	
geometry	str
radius	m
yheight	m
xwidth	m
zdepth	m
thickness	
pack	1
Vc	AA^3
sigma_abs	barns
sigma_inc	barns
delta_d_d	0/1
<b>p_inc</b>	1
<b>p_transmit</b>	1
<b>DW</b>	1
nb_scans	1
d_phi	deg
p_interact	
concentric	only for box, cylinder
density	g/cm^3
weight	g/mol
barns	1
Strain	ppm
focus_flip	1

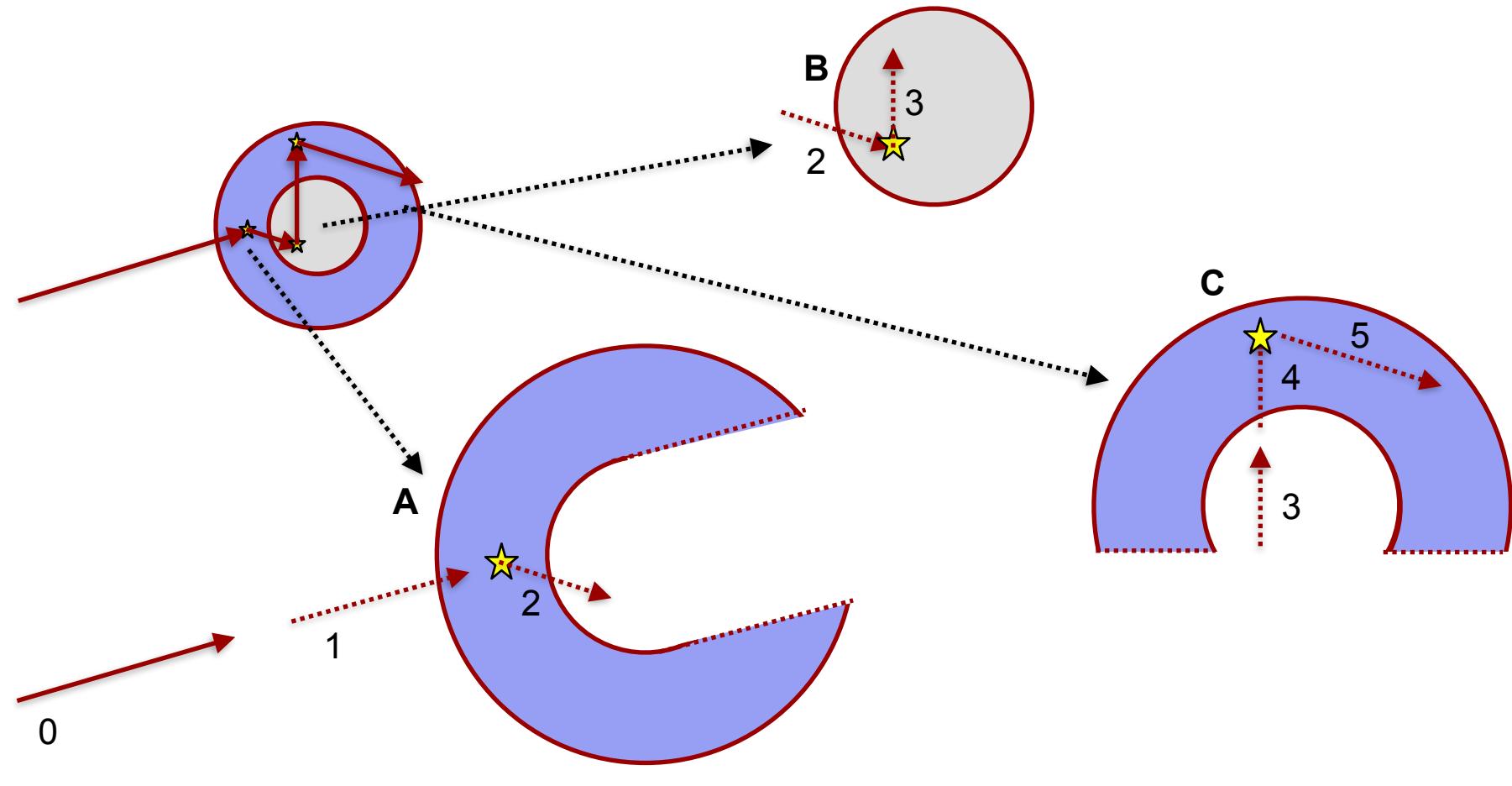**Concentric = Hollow  
Powder [cylinder]**


# PowderN concentric in code

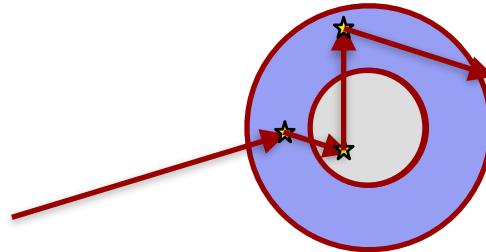

```
COMPONENT Al = PowderN(  
    reflections="Al.laz", radius=0.02, yheight=0.05, thickness=0.005,  
    p_transmit=0.8, concentric=1)  
AT (0, 0, 0) RELATIVE PREVIOUS
```

```
COMPONENT a12co3 = PowderN(  
    reflections=Al2CO3-laz, radius=0.005, yheight=0.05,  
    p_inc=0.15, p_transmit=0.5)  
AT (0, 0, 0) RELATIVE PREVIOUS
```

```
COMPONENT Al_2 = COPY(Al) (concentric=0)  
AT (0,0,0) RELATIVE PREVIOUS
```

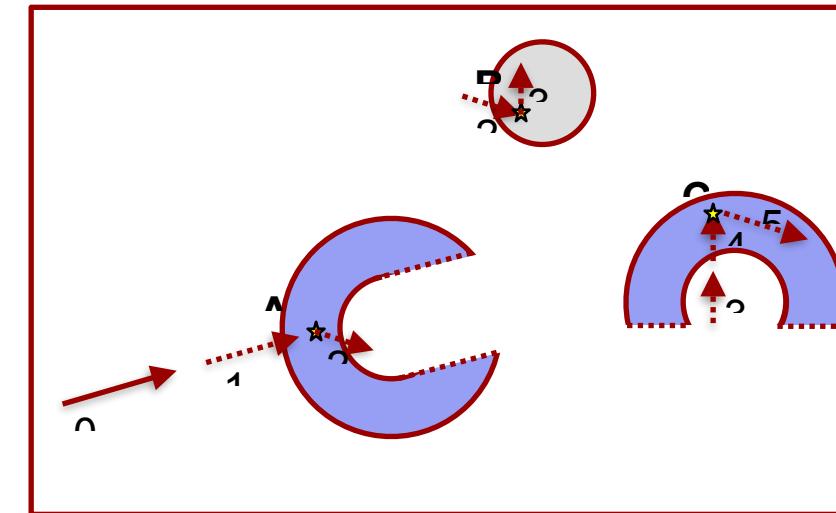

# PowderN inputs

## Input parameters

Parameters in **boldface** are required; the others are optional.

Name	Unit	Description	Default
format	no quotes	Name of the format, or list of column indexes (see Description).	Undefined
reflections			"NULL"
geometry	str	Name of an Object File Format (OFF) or PLY file for complex geometry. The OFF/PLY file may be generated from XYZ coordinates using qhull/powercrust	"NULL"
radius	m	Outer radius of sample in (x,z) plane	0
yheight	m	Height of sample y direction	0
xwidth	m	Horiz. dimension of sample, as a width	0
zdepth	m	Depth of box sample	0
thickness			0
pack	1	Packing factor	1
Vc	AA^3	Volume of unit cell=nb atoms per cell/density of atoms	0
sigma_abs	barns	Absorption cross section per unit cell at 2200 m/s. Use a negative value to unactivate it	0
sigma_inc	barns	Incoherent cross section per unit cell. Use a negative value to unactivate it	0
delta_d_d	0/1	Global relative delta_d_d/d broadening when the 'w' column is not available. Use 0 if ideal.	0
p_inc	1	Fraction of incoherently scattered neutron rays	0.1
p_transmit	1	Fraction of transmitted (only attenuated) neutron rays	0.1
DW	1	Global Debye-Waller factor when the 'DW' column is not available. Use 1 if included in F2	0
nb_atoms	1	Number of sub-unit per unit cell, that is ratio of sigma for chemical formula to sigma per unit cell	1
d_phi	deg	Angle corresponding to the vertical angular range to focus to, e.g. detector height. 0 for no focusing	0
p_interact	1	Fraction of events interacting with sample, e.g. 1-p_transmit-p_inc	0
concentric	only for box, cylinder, sphere	[1] Indicate that this component has a hollow geometry and may contain other components. It should then be duplicated after the inside part	0
density	g/cm^3	Density of material. rho=density/weight/1e24*N_A.	0
weight	g/mol		0
barns	1	s or fm^2 (barns=1 for laz, barns=0 for lau type files).	1
Strain	ppm	column is not available. Use 0 if ideal.	0
focus_flip	1	d against the xz-plane. If !=0 d_phi is measured against zy-plane.	0

Only scatter into a part of  
 $4\pi$


# PowderN inputs

## Input parameters

Parameters in **boldface** are required; the others are optional.

Name	Unit	Description	Default
format	no quotes	Name of the format, or list of column indexes (see Description).	Undefined
reflections			"NULL"
geometry	str	Name of an Object File Format (OFF) or PLY file for complex geometry. The OFF/PLY file may be generated from XYZ coordinates using qhull/powercrust	"NULL"
radius	m	Outer radius of sample in (x,z) plane	0
yheight	m	Height of sample y direction	0
xwidth	m	Horiz. dimension of sample, as a width	0
zdepth	m	Depth of box sample	0
thickness			0
pack	1	Packing factor	1
Vc	AA^3	Volume of unit cell=nb atoms per cell/density of atoms	0
sigma_abs	barns	Absorption cross section per unit cell at 2200 m/s. Use a negative value to unactivate it	0
sigma_inc	barns	Incoherent cross section per unit cell. Use a negative value to unactivate it	0
delta_d_d	0/1	Global relative delta_d_d/d broadening when the 'w' column is not available. Use 0 if ideal.	0
p_inc	1	Fraction of incoherently scattered neutron rays	0.1
p_transmit	1	Fraction of transmitted (only attenuated) neutron rays	0.1
DW	1	Global Debye-Waller factor when the 'DW' column is not available. Use 1 if included in F2	0
nb_atoms	1	Number of sub-unit per unit cell, that is ratio of sigma for chemical formula to sigma per unit cell	1
d_phi	deg	Angle corresponding to the vertical angular range to focus to, e.g. detector height. 0 for no focusing	0
p_interact	1	Fraction of events interacting with sample, e.g. 1-p_transmit-p_inc	0
concentric	only for box, cylinder, sphere	[1] Indicate that this component has a hollow geometry and may contain other components. It should then be duplicated after the inside part	0
density	g/cm^3	Density of material. rho=density/weight/1e24*N_A.	0
weight	g/mol		0
barns	1	s or fm^2 (barns=1 for laz, barns=0 for lau type files).	1
Strain	ppm	column is not available. Use 0 if ideal.	0
focus_flip	1	d against the xz-plane. If !=0 d_phi is measured against zy-plane.	0

Only scatter into a part of  
 $4\pi$

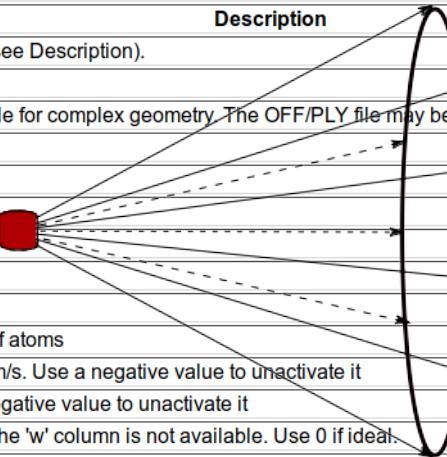

# PowderN inputs

## Input parameters

Parameters in **boldface** are required; the others are optional.

Name	Unit	Description	Default
format	no quotes	Name of the format, or list of column indexes (see Description).	Undefined
reflections			"NULL"
geometry	str	Name of an Object File Format (OFF) or PLY file for complex geometry. The OFF/PLY file may be generated from XYZ coordinates using qhull/powercrust	"NULL"
radius	m	Outer radius of sample in (x,z) plane	0
yheight	m	Height of sample y direction	0
xwidth	m	Horiz. dimension of sample, as a width	0
zdepth	m	Depth of box sample	0
thickness			0
pack	1	Packing factor	1
Vc	AA^3	Volume of unit cell=nb atoms per cell/density of atoms	0
sigma_abs	barns	Absorption cross section per unit cell at 2200 m/s. Use a negative value to unactivate it	0
sigma_inc	barns	Incoherent cross section per unit cell. Use a negative value to unactivate it	0
delta_d_d	0/1	Global relative delta_d_d/d broadening when the 'w' column is not available. Use 0 if ideal.	0
p_inc	1	Fraction of incoherently scattered neutron rays	0.1
p_transmit	1	Fraction of transmitted (only attenuated) neutron rays	0.1
DW	1	Global Debye-Waller factor when the 'DW' column is not available. Use 1 if included in F2	0
nb_atoms	1	Number of sub-unit per unit cell, that is ratio of sigma for chemical formula to sigma per unit cell	1
d_phi	deg	Angle corresponding to the vertical angular range to focus to, e.g. detector height. 0 for no focusing	0
p_interact	1	Fraction of events interacting with sample, e.g. 1-p_transmit-p_inc	0
concentric	only for box, cylinder, sphere	[1] Indicate that this component has a hollow geometry and may contain other components. It should then be duplicated after the inside part	0
density	g/cm^3	Density of material. rho=density/weight/1e24*N_A.	0
weight	g/mol		0
barns	1	s or fm^2 (barns=1 for laz, barns=0 for lau type files).	1
Strain	ppm	column is not available. Use 0 if ideal.	0
focus_flip	1	d against the xz-plane. If !=0 d_phi is measured against zy-plane.	0

Only scatter into a part of  
 $4\pi$

# Reflection files for Single\_crystal and PowderN

```
# TITLE *Aluminum-Al-[FM3-M] Miller, H.P.jr.;DuMond, J.W.M.[1942] at 298 K
# CELL 4.049320 4.049320 4.049320 90.000000 90.000000 90.000000
# SPCGRP F M 3 M CUBIC STRUCTURE
# ATOM AL 1 0.000000 0.000000 0.000000
# SCATTERING FACTOR COEFFICIENTS: AL F= 0.345 CM-12
# Reference: Physical Review (1940) 57, 198-206
#
# Physical parameters:
# sigma_coh 1.495 coherent scattering cross section (single atom) in [barn]
# sigma_inc 0.0082 incoherent scattering cross section (single atom)in [barn]
# sigma_abs 0.231 absorption scattering cross section (single atom) in [barn]
# density 2.70 in [g/cm^3]
# weight 26.98 in [g/mol] (single atom)
# multiplicity 4 in [atoms/unit cell]
# Vc 66.4 volume of unit cell in [A^3]
# v_sound 5100 in [m/s]
# v_sound_l 6420 velocity of longitudinal sound in [m/s]
# v_sound_t 3040 velocity of transversal sound in [m/s]
# T_m 933.5 melting temperature in [K]
# T_b 2792.2 boiling temperature in [K]
# At_number 13 atomic number Z
# lattice_a 4.04932 lattice parameter a in [Angs]
#
# Format parameters: Crystallographica format
# column_j 4 multiplicity 'j'
# column_d 5 d-spacing 'd' in [Angs]
# column_F2 7 norm of scattering factor |F|^2 in [fm^2]
# column_h 1
# column_k 2
# column_l 3
#
# h k l Mult. d-space 2Theta F-squared
-1 -1 8 2.338 24.6973 21.3
-1 -1 1 8 2.338 24.6973 21.3
...

```

Lau datafiles

header

+

reflection list

Can be used with  
Single\_crystal, PowderN,  
Isotropic\_Sqw

# Input for the PowderN

```

# TITLE *Corundum-Al2O3-[R3-CH] Graafsma, H.;Souhassou, M.;Harkem[1998] [corundum saphire:blue, ruby:red]
# CELL 4.757000 4.757000 12.987700 90.000000 90.000000 120.000000
# SPCGRP R -3 C TRIGONAL STRUCTURE
# ATOM AL 1 0.000000 0.000000 0.352110
# ATOM O 1 0.306260 0.306260 0.250000
# SCATTERING FACTOR COEFFICIENTS: AL F= 0.345 CM-12 ; O F= 0.581 CM-12
# Reference: Acta Crystallographica B (1998) 54, 193-195
#
# Physical parameters:
# sigma_coh 15.683 coherent scattering cross section for Al2O3 in [barn]
# sigma_inc 0.0188 incoherent scattering cross section for Al2O3 in [barn]
# sigma_abs 0.4625 absorption scattering cross section for Al2O3 in [barn]
# density 4.05 in [g/cm^3]
# weight 101.96 in [g/mol] for Al2O3
# multiplicity 6 in [Al2O3/unit cell]
# Vc 254.52 volume of unit cell in [A^3]
# T_m 2273 melting temperature in [K]
# T_b 3773 boiling temperature in [K]
# lattice_a 4.757 lattice parameter a in [Angs]
# lattice_c 12.9877 lattice parameter c in [Angs]
# lattice_cc 120 lattice angle gamma in [deg]
#
# Format parameters: Lazy format <http://icsd.ill.fr>
# column_j 17 multiplicity 'j'
# column_d 6 d-spacing 'd' in [Angs]
# column_F 13 norm of scattering factor |F| in [barn]
# column_h 1
# column_k 2
# column_l 3
#
# H K L THETA 2THETA D VALUE 1/D**2 SIN2*1000 H K L INTENSITY /F(HKL)/ A(HKL) B(HKL) PHA.ANG. MULT LPG
1 0 1 6.35 12.71 4.5175 0.0490 12.25 1 0 1 367.0 4.1 -4.08 0.00 180.00 6 82.14
0 0 3 7.10 14.20 4.0467 0.0611 15.27 0 0 3 110.0 4.3 4.32 0.00 0.00 2 66.01
0 1 2 7.57 15.13 3.7972 0.0694 17.34 0 1 2 10.9 0.8 0.84 0.00 0.00 6 58.18
...

```

Laz + Lau datafiles

header

+

reflection list

Can be used with PowderN,  
Isotropic\_Sqw

# What's already there?

## Single Crystals

Al.lau  
 Al203\_sapphire.lau  
 B4C.lau  
 BCC\_fake.lau  
 BeO.lau  
 C60.lau  
 C\_diamond.lau  
 C\_graphite.lau  
 CaF2.lau  
 Ge.lau  
 In.lau  
 PPase\_D\_P1.lau  
 Rubredoxin.lau  
 SiO2\_cristobalite.lau  
 SiO2\_quartz.a.lau  
 SiO2\_quartzb.lau  
 TiO2\_rutile.lau  
 YBaCuO.lau  
 adrenaline.lau  
 aspirin.lau  
 leucine.lau


## Powders

Ag.laz  
 Al.laz  
 Al203\_sapphire.laz  
 Au.laz  
 B4C.laz  
 Ba.laz  
 Be.laz  
 BeO.laz  
 Bi.laz  
 C\_diamond.laz  
 C\_graphite.laz  
 Cd.laz  
 Cr.laz  
 Cs.laz  
 Cu.laz  
 Cu2MnAl.laz  
 Fe.laz  
 Fe\_Gamma.laz  
 Ga.laz  
 Gd.laz  
 Ge.laz  
 H2O\_ice\_1h.laz  
 He4\_hcp.laz  
 Hg.laz  
 I2.laz  
 K.laz  
 Li.laz  
 LiF.laz  
 Mo.laz  
 Na2Ca3Al2F14.laz  
 Nb.laz  
 Ni.laz  
 Pb.laz  
 Pd.laz  
 Pt.laz  
 Rb.laz  
 Se\_alpha.laz  
 Se\_beta.laz  
 Si.laz  
 SiO2\_quartz.a.laz  
 SiO2\_quartzb.laz  
 Sn\_alpha.laz  
 Sn\_beta.laz  
 Ti.laz  
 Tl.laz  
 UO2.laz  
 V.laz  
 Y2O3.laz  
 Y3Fe5O12\_YIG.laz  
 Zn.laz  
 Zr.laz  
 duplex.laz

# Where to get these files...

- | \$MCSTAS/data
  - | Windows: c:\mcstas-3.1\lib\data
  - | Linux: /usr/(local)/share/mcstas/3.1/data
  - | OS X: /Applications/McStas-3.1/Contents/Resources/mcstas/3.1/data
- | - Or make your own via
  - | Finding a CIF file for the given structure
    - | e.g. from  
ICSD <http://icsd.fiz-karlsruhe.de> or  
COD <<http://crystallography.net>>
- | Process it using
  - | cif2hkl which is available in your McStas installation

[https://github.com/McStasMcXtrace/McCode/wiki/Generating-input-files-for-Single\\_crystal-and-PowderN](https://github.com/McStasMcXtrace/McCode/wiki/Generating-input-files-for-Single_crystal-and-PowderN)

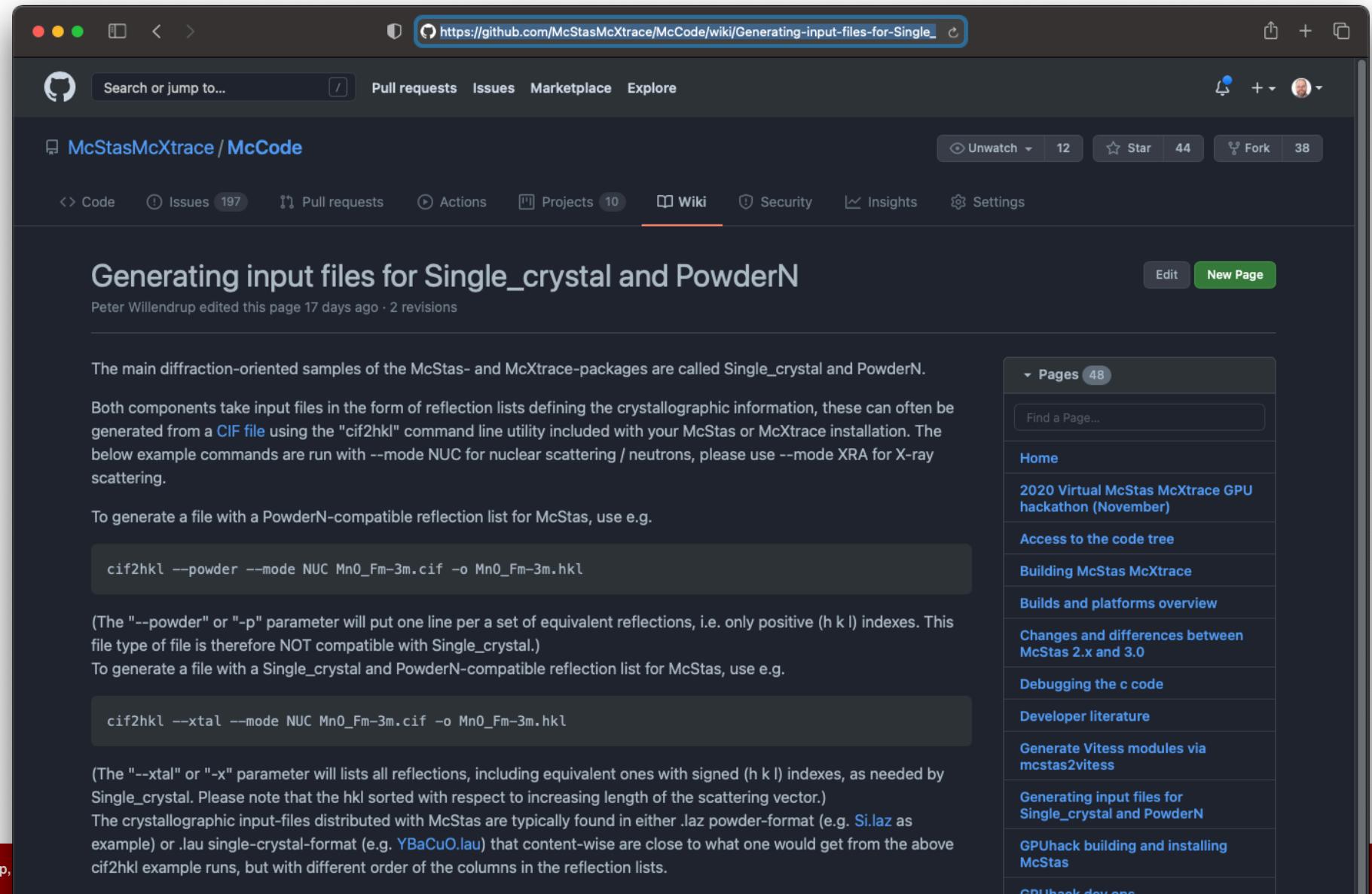

The screenshot shows a GitHub repository page for 'McStasMcXtrace / McCode'. The specific page displayed is 'Generating input files for Single\_crystal and PowderN'. The page content provides instructions on how to generate reflection lists for McStas and McXtrace. It includes examples of command-line usage for 'cif2hkl' and 'cif2hkl --xtal'. A sidebar on the right lists other pages in the wiki, such as 'Access to the code tree', 'Building McStas McXtrace', and 'Changes and differences between McStas 2.x and 3.0'.

**Generating input files for Single\_crystal and PowderN**

Peter Willendrup edited this page 17 days ago · 2 revisions

The main diffraction-oriented samples of the McStas- and McXtrace-packages are called Single\_crystal and PowderN.

Both components take input files in the form of reflection lists defining the crystallographic information, these can often be generated from a CIF file using the "cif2hkl" command line utility included with your McStas or McXtrace installation. The below example commands are run with --mode NUC for nuclear scattering / neutrons, please use --mode XRA for X-ray scattering.

To generate a file with a PowderN-compatible reflection list for McStas, use e.g.

```
cif2hkl --powder --mode NUC Mn0_Fm-3m.cif -o Mn0_Fm-3m.hkl
```

(The "--powder" or "-p" parameter will put one line per a set of equivalent reflections, i.e. only positive (h k l) indexes. This file type of file is therefore NOT compatible with Single\_crystal.)

To generate a file with a Single\_crystal and PowderN-compatible reflection list for McStas, use e.g.

```
cif2hkl --xtal --mode NUC Mn0_Fm-3m.cif -o Mn0_Fm-3m.hkl
```

(The "--xtal" or "-x" parameter will lists all reflections, including equivalent ones with signed (h k l) indexes, as needed by Single\_crystal. Please note that the hkl sorted with respect to increasing length of the scattering vector.)



The crystallographic input-files distributed with McStas are typically found in either .laz powder-format (e.g. Si.laz as example) or .lau single-crystal-format (e.g. YBaCuO.lau) that content-wise are close to what one would get from the above cif2hkl example runs, but with different order of the columns in the reflection lists.

Pages 48

- Find a Page...
- Home
- 2020 Virtual McStas McXtrace GPU hackathon (November)
- Access to the code tree
- Building McStas McXtrace
- Builds and platforms overview
- Changes and differences between McStas 2.x and 3.0
- Debugging the c code
- Developer literature
- Generate Vitess modules via mcstas2vitess
- Generating input files for Single\_crystal and PowderN
- GPUhack building and installing McStas
- GPUhack dev doc

# McStas: structure: powder

- For powders, the pure elastic diffraction is easy to prepare: LAZ
  - PowderN
  - Single\_crystal(powder=1, ...)
  - Isotropic\_Sqw

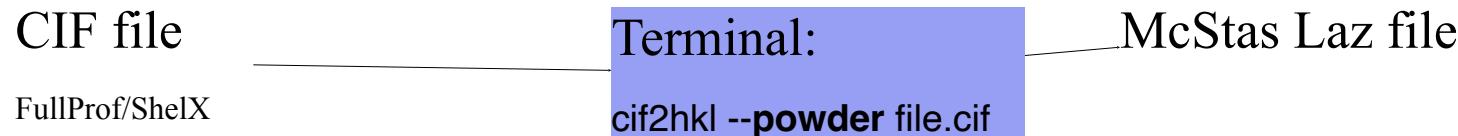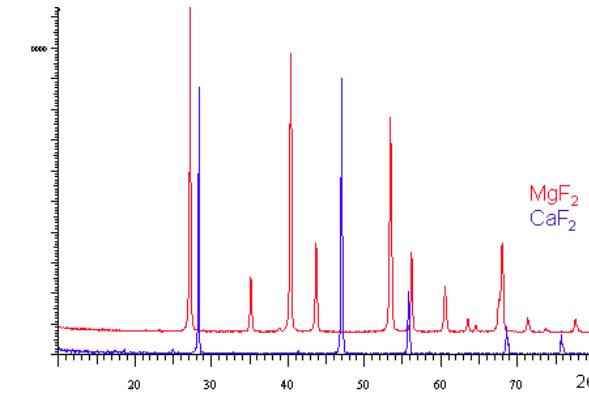

# McStas: structure: SX

- For single crystals, the pure elastic diffraction is easy to prepare: LAU
  - PowderN
  - Single\_crystal
  - Isotropic\_Sqw



CIF file

FullProf/ShelX

Terminal:

```
cif2hkl --xtal file.cif
```

McStas Lau file

## McStas Data formats

McStas uses a number of input data formats, all text based.

- **Lau** (Single\_crystal, PowderN, Isotropic\_Sqw)
- **Laz** (PowderN, Isotropic\_Sqw)
- **Sqw** (Isotropic\_Sqw)
- **qSq** (Isotropic\_Sqw) – mostly for liquids
- **PDB** (SANS\_PDB)
- **ncmat** (NCrystal)

# OFF/PLY any geometry

An additional complex geometry enables to use any point set to describe the material volume (geomview OFF file).

